

Stem Cell Exosome Production on the SciVario® twin, a Flexible Controller for Your Bioprocess Needs.

Jorge L. Escobar Ivirico and Ma Sha

Eppendorf, Inc., Enfield, CT, USA

Contact: bioprocess-experts@eppendorf.com

Abstract

Exosomes derived from mesenchymal stem cells (MSC) are a new alternative to restore tissues and organs, avoiding the limitations associated with stem cell therapy. MSC-derived exosomes are involved in MSC's paracrine functions related to cell-to-cell communications and tissue reconstruction, leading to cell self-renewal and differentiation. Despite the benefits of stem cell use, there are limitations, including the large-scale production of cells that mimic physiological conditions. Therefore, increasing the mass production of exosomes in a controlled environment is the necessary next step. In this study, we used a combination of SciVario® twin bioreactor control system and BioBLU® 1c Single-Use Bioreactors as a model for large scale exosome production in bioreactors. The

SciVario twin is a future-proof bioreactor control system developed by Eppendorf. The flexibility of the controller allows working volumes between 1.3 and 40 L, depending on the bioreactor type. Following the agile development principle, further hard- and software updates are possible making the SciVario twin a good choice for current and future bioprocessing needs.

In this application note, we monitored and analyzed the metabolites derived from the human induced Pluripotent Stem Cells (hiPSC) as well as cell adhesion and morphology, suggesting that we obtained an excellent physiological environment over time. Furthermore, we isolated, purified and quantified (by CD63 ELISA kit) the exosomes produced from MSC.

Introduction

Regenerative medicine is a multidisciplinary field that engineers the structure and function of tissues and organs. Due to their ability to migrate to the site of injury and promote tissue regeneration through paracrine factors (secretome), mesenchymal stem cells have become the most widely used stem cell type for such investigations [1-3]. However, problems associated with inadequate cell localization and low cell survival rate within the target tissue make the MSC less attractive. Recently, paracrine factors have attracted increasing interest due to their potential at overcoming the limitations of MSC. Extracellular vesicles (EV), including exosomes, are one of the most important paracrine effectors involved in intracellular communication and trafficking [4].

Exosomes are lipid bilayer vesicles with a diameter ranging from 30 to 200 nm, which can be easily confirmed by surface

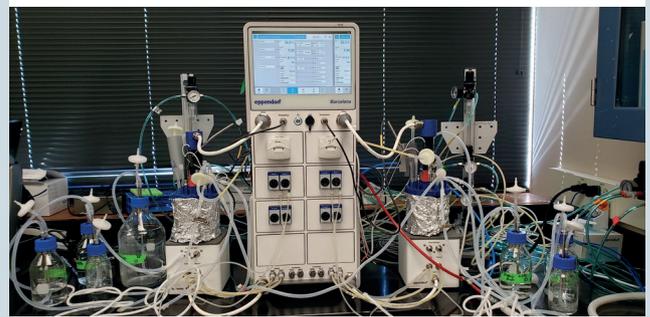


Fig. 1: Human iPSC-derived MSC exosome production in a medium exchange-batch culture using the SciVario twin bioreactor control system and BioBLU 1c Single-Use Bioreactors.

markers such as CD9, CD63, CD81 of tetraspanins [5]. They offer several advantages compared to cell-based approaches, including exceptional stability and biocompatibility. Moreover, they can be delivered to target tissues easier than whole cells, and can migrate across the blood brain barrier [6]. In addition, since exosomes have the ability to conduct immune modulation, there is much less risk for tumor formation or inhibition of inflammatory cell migration [7, 8]. Moreover, the lack of complex metabolism of exosomes and the influence of the environment on target tissues reduces complications upon use.

Material and Methods

SciVario twin

We used the SciVario twin bioreactor control system to perform medium exchange-batch culture using BioBLU 1c Single-Use Bioreactors equipped with a single pitched-blade impeller. Each bioreactor unit possesses three universal port connectors for pH (port 1) and DO (dissolved oxygen; port 2) sensors, a temperature control block that combines electrical heating and water cooling, agitation control, and a gas module that includes 1 TMFC (Thermal Mass Flow Controller) with an ultra-high turndown ratio of 1:12,000, and 4 solenoid valves for automated 4-gas mixing.

Cell culture in T-flasks and multilayers flasks

We purchased the iPSC-derived MSC from American Type Culture Collection (ATCC®, ACS-7010) and we established the culture by seeding 10,000 cells/cm² per T-75 flask using either 15 mL of DMEM/F12 medium (Thermo Fisher Scientific®, 11320033) supplemented with 10 % fetal bovine serum (FBS) (Thermo Fisher Scientific, 16000044) and 1 % antibiotic-antimycotic (100X) (Sigma Aldrich, A5955) or mesenchymal stem cell basal medium (ATCC, PCS-500-030) (ATCC medium) supplemented with 7 % FBS, (Thermo Fisher Scientific, 16000044), 125 pg/mL rhFGF basic, 15 ng/mL rh IGF-1, 2.4 mM L-alanyl-L-glutamine (Mesenchymal Stem Cell Growth Kit (ATCC PCS-500-041) and 1 % antibiotic-antimycotic (100X) (complete ATCC medium).

We incubated the flasks at 37 °C and 5 % CO₂ in a CellXpert® C170i Incubator. When the cultures had reached approximately 80 % to 90 % confluence, we made a second passage as follow: For each flask, we aspirated the spent me-

di-um, rinsed the cell layer with 5 mL of Dulbecco's phosphate buffered saline (DPBS) (Thermo Fisher Scientific, 14190144) to remove residual traces of serum and aspirated the DPBS. Then, we added 3 mL of pre-warmed trypsin-EDTA solution (0.25 %) (Thermo Fisher Scientific, 25200056) and incubated flasks at 37 °C and 5 % CO₂ for 5 minutes. When most cells were detached, we added the double volume of the ATCC complete medium to each flask. We collected and counted the cells using a Vi-CELL XR cell viability analyzer (Beckman Coulter®). Finally, we centrifuged the iPSC-derived MSCs at 1500 rpm (Centrifuge 5430R) for 5 minutes (we aspirated the neutralized dissociation solution from the cell pellet), resuspended in DMEM/F12 medium or ATCC complete medium and seeded at a density of 10,000 cells/cm² in T-175 flasks.

We performed a third passage (required to achieve the appropriate inoculation cell density) following the same cell subculture procedure described above, but seeding the cells in multilayer flasks, a HYPERFlask® M cell culture vessel (Corning®, 09-761-22) with the same cell culture density.

We performed a third passage (required to achieve the appropriate inoculation cell density) following the same cell subculture procedure described above, but seeding the cells in multilayer flasks, a HYPERFlask® M cell culture vessel (Corning®, 09-761-22) with the same cell culture density.

Microcarrier preparation

We used collagen-coated cross-linked polystyrene microcarriers (Pall Corporation, C-221-020) as support matrix. We conducted the sterilization process according the manufacturer's instructions [10] but with a few modifications. Briefly, we transferred 17 g of collagen coated microcarriers (125 – 212 µm of size range and 360 cm²/g of superficial area) into a 250 mL glass bottle along with 100 mL of DPBS. Then, we autoclaved the glass bottle at 121°C for 30 minutes.

After the sterilization process, we carefully aspirated the autoclaved PBS, added 100 mL of medium to the microcarriers and incubated for 24 hours.

Sensor calibration

Prior to the preparation of the BioBLU 1c Single-Use Bioreactors, we connected the ISM® gel-filled pH sensors (Mettler Toledo®) to the SciVario twin bioreactor control system that was automatically detected by the software of the controller. We performed the calibration process according to the operation’s manual using buffer solutions of pH 7 and pH 4 as “zero” and “span” respectively. Then, we disconnected the pH sensors and sterilized them in an autoclavable pouch.

BioBLU 1c Single-Use Bioreactor preparation and process parameters

Each BioBLU 1c Single-Use Bioreactors was equipped with a pH sensor, as well as a dip tube along with a compression fitting adapter, both inserted in a spare Pg 13.5 port under aseptic conditions in the BioSafety Cabinet following sterilization. In addition, we fitted a polarographic DO sensor (Mettler Toledo) in the headplate, an exhaust condenser, a sampling dip tube, a 4-gas mixing line connected to the gas sparge port and an extra N₂ gas line connected to an overlay gas port and 3 liquid addition ports. There is one additional port for the inoculation and medium addition, one for the addition of the base and another for the addition of 0.1 % antifoam (Pluronic®-F68 surfactant, Life Technologies®, 24040-032). We then placed the assembled vessels in their respective temperature control blocks to maintain constant temperature. Finally, we introduced the DMEM/F12 or ATCC complete medium into each vessel and conditioned the vessels for at least 24 hours under the parameters and setpoints listed in Table 1.

Cascade control of DO

We establish the following DO cascade to control the culture with air only, without oxygen supplementation:
 “Set O₂ % at 30 % controller output to 21 % and at 100 % controller output to 21 %. Set flow at 0 % controller output to 0.50 SLPH and at 100 % controller output to 30 SLPH.”

hiPSC-derived MSCs culture on BioBLU 1c Single-Use Bioreactors

We harvested the cells from multilayer flasks (see section “Cell culture in T-flasks and multilayers flasks”) and seeded at a density described in Table 1 into the glass bottles containing collagen-coated microcarriers (see section “Microcarriers preparation”) and 300 mL of DMEM/F12 or ATCC complete medium. We kept the bottles in a CellXpert C170i

Table 1: Process parameters and setpoints of the first and second experiments.

	First Experiment	Second Experiment
Parameters	Setpoints	
Starting volume	700 mL	
Ending volume	1 L	
Initial agitation	80 rpm (0.2 tip speed)	
Temperature	37 °C	
Inoculation density	3 x 10 ⁴ cell/mL	10.4 x 10 ⁴ cells/ mL
Cell culture medium	DMEM/F12 medium	ATCC complete medium
DO Setpoint	40% (P=0.1; I=3.6/h)	
pH Setpoint	7.2 (deadband = 0.1), cascade to CO ₂ (acid) cascade to 0.45 M sodium bicarbonate (base)	7.6 (deadband = 0.1), cascade to CO ₂ (acid) cascade to 0.45 M sodium bicarbonate (base)
Overlay N₂ gas flow	0.20 SLPM	0.25 SLPM
Gassing range	0.1 SLPH-30 SLPH	
Gassing cascade	Set O ₂ % at 30 % controller output to 21 % and at 100 % controller output to 21 %. Set flow at 0 % controller output to 0.1 SLPH, and at 100 % controller output to 30 SLPH.	

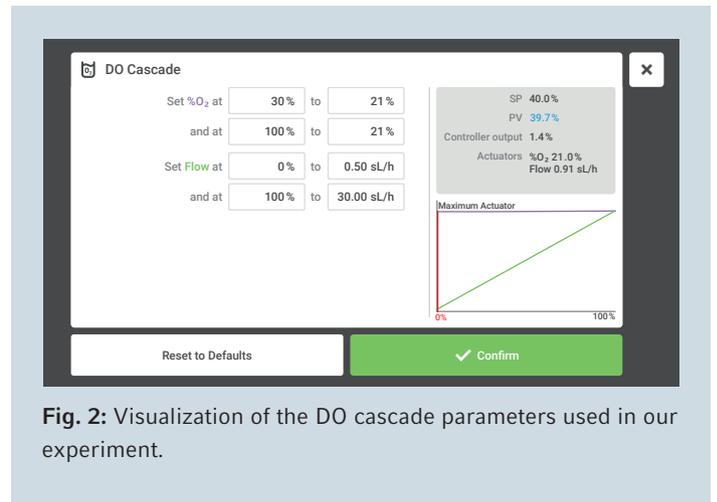


Fig. 2: Visualization of the DO cascade parameters used in our experiment.

incubator for 2 hours at 37 °C and 5 % CO₂ under static conditions. Following incubation, the 300 mL of cells on microcarriers were then transferred to a BioBLU 1c Single-Use Bioreactor, already filled with 700 mL medium, reaching the working volume of 1 L and maintaining the microcarriers concentration at 17 g/L. Furthermore, we set the overlay N₂ gas flow between 0.2 and 0.25 SLPM to maintain the DO level at 40 %.

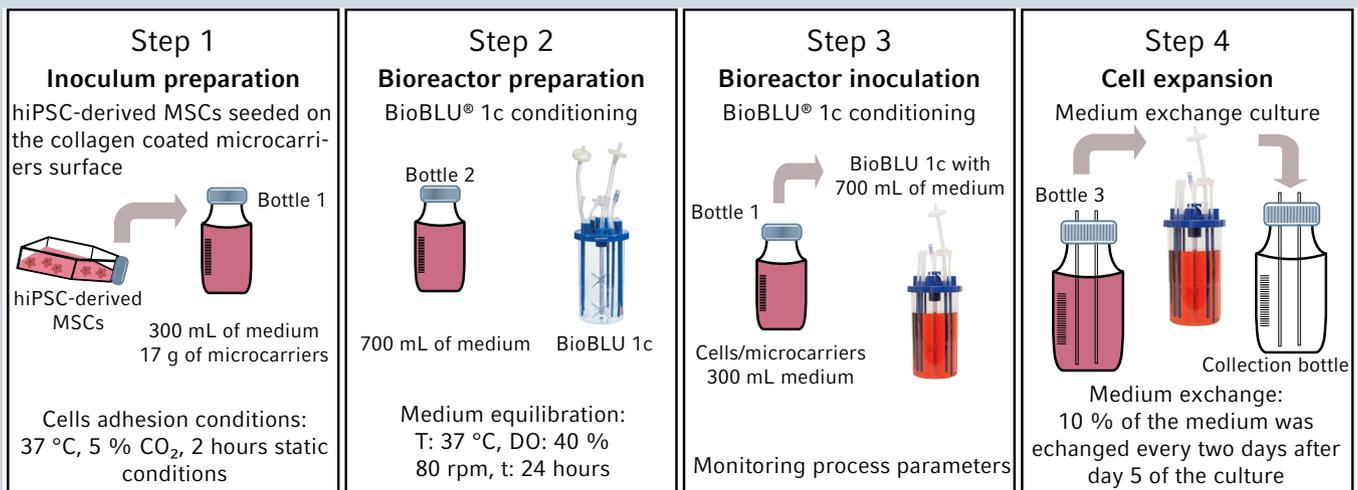


Fig 3: Process flow diagram of the experimental setup

Medium exchange

As we previously described in the “**BioBLU 1c Single-Use Bioreactor preparation and process parameters**” section, we inserted a dip tube along with a compression probe adapter in a spare Pg 13.5 port allowing for medium exchange without disturbing the cell culture. Briefly, after day 5 of cell culture, we performed a medium exchanged every two days until day 9 and then on daily basis as required. We stopped the agitation and gases flow for 5 minutes, allowing the microcarriers to settle at the bottom of the BioBLU 1c and exchanged 10 % of medium, using the dip tube to remove the cell culture medium from the surface and using the feeding port to add fresh medium. Moreover, in the second experiment, additional glucose solution with a concentration of 1 g/L was added on day 11 of the culture to the fresh medium to increase glucose levels in the cell culture with subsequent medium exchanges.

Exosome collection assay


After day 5 of culture, we collected 50 mL of iPCS-derived MSC/microcarriers/medium using a Labtainer Pro BioProcess Container (BPC) bag with line sets (Thermo Fisher Scientific, SH30658.11). We transferred the 50 mL into 50 mL conical tube allowing the iPCS-derived MSC/microcarriers system to settle to the bottom of the tube. We removed the supernatant and added additional 50 mL of medium with the same composition (dependent on the experiment) but replacing the regular FBS by 2 % exosome-depleted FBS (ED-FBS) (Fisher Scientific, A2720801) and transferring the contents into a 250 mL shake flask (Schott Duran®). Finally, we incubated

the iPSC-derived MSC/microcarriers system for 48 hours at 37 °C, 5 % CO₂ and 50 rpm using the New Brunswick S41i CO₂ Incubator Shaker. We repeated the collection process at days 8, 11, and 14 of culture.

Isolation and purification of human iPSC-derived MSC exosomes

At 48 hours after each collection day, we removed the supernatant and performed the isolation and purification of exosomes according to the ExoQuick-TC PLUS (System Bioscience, EQPL10TC-1) protocol with some modifications. Briefly, we centrifuged the supernatant at 3000 × g (Centrifuge 5430R) for 15 minutes and filtered using a 0.22 μm filter to remove cells and cell debris. Then, we filtered the supernatant using the Amicon® Ultra-15 Centrifugal system (Millipore Sigma, UFC910024), with molecular weight cut-off (MWCO) of 100,000, to concentrate the exosome fraction. We transferred the supernatant to a sterile tube, added ExoQuick-TC Exosome Precipitation Solution in a ratio of 5:1 and incubated at 4 °C overnight. After the incubation period, we centrifuged the ExoQuick-TC/supernatant mixture at 1500 × g for 30 minutes, bringing the exosomes pellet to the bottom of the tube. Then, we carefully aspirated the supernatant, centrifuged the pellet again at 1500 × g for 5 minutes to remove any residue of ExoQuick-TC solution and resuspended the exosomes in 250 μL of resuspension buffer.

Finally, we washed the microsphere beads three times with resuspension buffer and added 250 μL of the exosome sample to the washed microspheres. The microspheres beads are customized to reduce protein carryover from serum/plasma

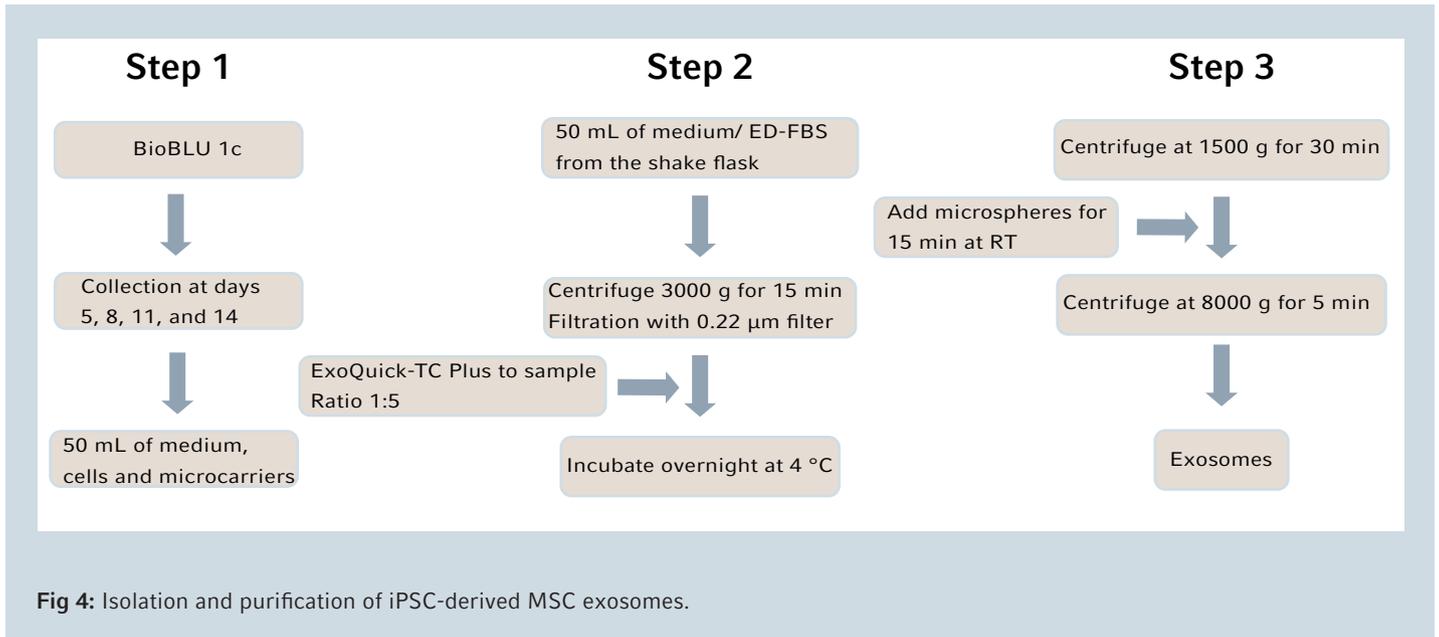

and tissue culture. We collected the supernatant containing the purified exosomes after mixing the exosomes and microspheres at room temperature on an inverting shaker for 15 minutes and centrifuged the mixture for 5 minutes at 8500 rpm (Centrifuge MiniSpin® Plus).

Cell viability and metabolic activity

We collected a sample every day in the course of the first experiment or every two days in the second experiment from the bioreactors to determine the cell viability, cellular density, and metabolites concentration (glucose, ammonia (NH₃), and lactate). This was accomplished by connecting a sterile 5 ml syringe to the Luer Lock sample port. Then, we discarded 3 mL of dead volume and collected again 3 mL (using a new 5 mL sterile syringe) as a workable sample for analysis. We transferred the viable samples to a 15 mL Eppendorf conical tube, allowing the microcarriers to settle in the tube. Before aspirating the medium, we used 1 mL to measure the metabolite levels employing a Cedex® Bio Analyzer (Roche). Then, we washed the microcarriers with 1 mL of DPBS and incubated at 37 °C for 5 minutes with 300 µL of trypsin (0.25%). To inactivate the trypsin solution, we added 700 µL of fresh medium and gently pipetted the sample up and down, creating a cell suspension. We filtered the cells and microcarriers

through a 40-micron strainer into a 50 mL Eppendorf conical tube. Finally, we measured the cell viability and density using a Vi-CELL XR Viability Analyzer.

Quantification of human iPSC-derived MSC exosomes

Following the exosomes enrichment by the ExoQuick-TC PLUS protocol, we quantified tetraspanin-containing exosomes through ExoELISA-ULTRA CD63 Kit (System Bioscience, EXEL-ULTRA-CD63-1). We immobilized the exosomes on the surface of the CD63 ExoELISA microtiter plate. Then, we finished the assay according to the manufacture's protocol instructions.

iPSC-derived MSC morphology

To study the human iPSC-derived MSC morphology during the cell expansion process in the BioBLU 1c Single-Use Bioreactor, we followed the protocol described below: We collected 500 µL of cells/microcarriers/medium at days 5 and 9 of culture and transferred them to a 1.5 mL Eppendorf microtube, allowing the microcarriers to settle to the bottom. Then, we removed most of the supernatant and used the cells and microcarriers to facilitate visualization of the bright-field images, using an EVOS® FL fluorescence microscope (Life Technology).

Results and Discussion

To demonstrate the ability of human iPSC-derived MSC to secrete exosomes, we performed two medium exchanges using the BioBLU 1c Single-Use Bioreactors and SciVario twin bioreactor control system. After the initial expansion of the cells in T-flasks culture conditions, we analyzed their stemness capacity by flow cytometry noting that they displayed the typical phenotype of MSCs (see Figure 5). Specifically, cells were positive for CD90 and CD29 (typical MSC markers) and negative for hematopoietic markers such as CD34 and CD11b.

Next, we established the optimal cell culture conditions. We used 17 g/L microcarrier concentration to provide cell support and employing DMEM/F12 as the cell culture medium. The inoculum was ready after the cell expansion and the initial cell attachment process provided by the microcarrier in the CellXpert C170i incubator (at 37 °C and 5 % CO₂ for 2 hours under static conditions). We then inoculated the BioBLU 1c Single-Use Bioreactor at an initial cell density of 5 × 10³ cells/cm² (3 × 10⁴ cell/mL) under a controlled environment. In addition, we added the anti-foaming agent Pluronic-F68 surfactant (0.1 %) to the medium.

After the initial cell expansion, the cell count was erratic, especially after each collection day (Figure 6). We attributed this behavior to the low inoculation density and the use of

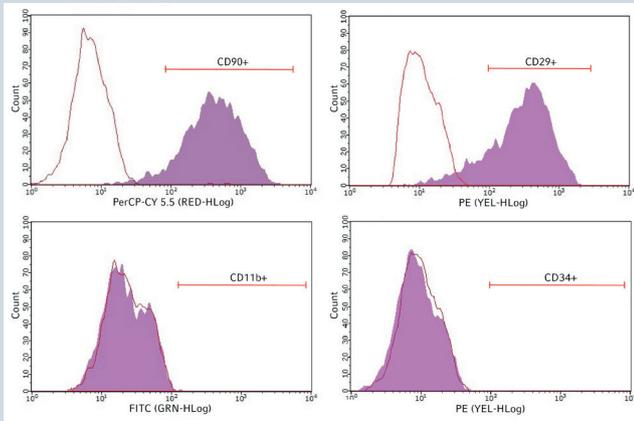


Fig. 5: Immunophenotyping of human iPSC-derived MSC at passage 3. We immunophenotyped the human iPSC-derived MSC for CD90, CD29, CD34, and CD11b (filled histogram). We used unstained cells as control (empty histogram).

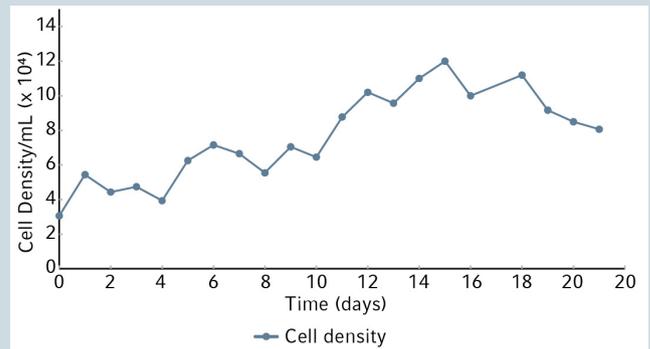


Fig. 6: iPSCs-derived MSCs cell growth trending BioBLU 1c Single-Use Bioreactor with DMEM/F12 medium.

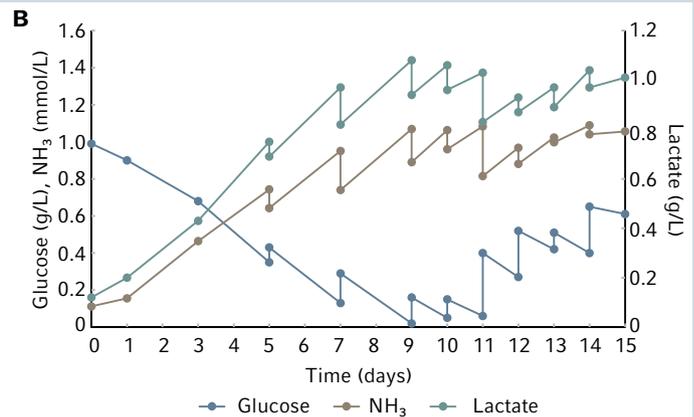
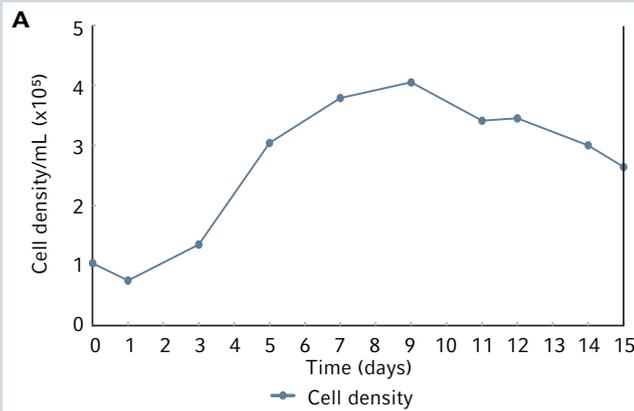


Fig. 7: iPSCs-derived MSCs growth profile in BioBLU 1c Single-Use Bioreactor with ATCC complete medium. **(A)** iPSCs-derived MSC density and viability. **(B)** Metabolic profile.

DMEM/F12 cell culture medium whose composition may be suboptimal for the iPSC-derived MSC expansion. Additionally, the DO level is known to be an important factor that can affect the growth and performance of the MSC. Although some studies reported that hypoxic conditions enhance the therapeutic potential of MSC and its secretome, others indicated that no matter under what conditions the experiment is established (hypoxia or normoxia), the collected secretome showed similar effects [11]. In this study, the O₂ demand is expected to be low due to the low inoculation cell density used, which is why we established a cascade with air flow only, without oxygen supplementation (21%) (Figure 2) throughout the experiment to reach the desired setpoint. We also added an N₂ overlay to balance the system using an external flow meter (Omega Engineering, FMA-1608A) directly connected to the nitrogen gas outlet, bypassing the controller. Overall, cell growth increased 4-fold more than the initial cell density on day 15 of culture, but the final density was low.

We performed a second experiment to increase the iPSCs-derived MSCs inoculation density and the exosomes production yield from each sample. To achieve that, we changed the cell culture conditions of the iPSC-derived MSC relative to the first experiment. First, we performed the cell expansion and the initial cell attachment to the microcarriers in ATCC complete medium. Then, we increased the inoculation cell density to 17×10^3 cells/cm² (10.4×10^4 cell/mL), maintaining the microcarrier concentration at 17 g/L. We kept the rest of the conditions the same as for the first experiment, with the exception that the overlay gas flow was increased to 0.25 SLPM to provide a better balance to the minimum air flow. Additionally, we performed a medium exchange on a daily basis starting from day 5 of the culture. We observed an initial lag phase 24 hours after the inoculation followed by a steady increase of cell growth between days 1 and 9 of culture.

We then determined the consumption of glucose and production of lactate and NH₃ while maintaining the concentration of lactate and NH₃ below 1.2 g/L and 1.2 mmol/L respectively during the whole run. The glucose level was significantly lower in the ATCC complete medium than in the DMEM/F12, necessitating the addition of 1 g/L of glucose.

The culture reached the stationary phase around day 9, followed by a subsequent decrease in cell density until day 15 (Figure 7). At day 9, the iPSC-derived MSC reached a maximum cell density of 4.1×10^5 cell/mL. However, we observed substantial microcarrier aggregation later in the course of the run. Therefore we based the cell count in the later stages on only the floating microcarriers, thus the true average of the cell expansion in the vessel is undoubtedly higher.

To isolate and characterize the iPSC-derived MSC exosomes, we used a labtainer bag with line sets connected to the BioBLU 1c. At day 5, 8, 11 and 14, we collected the exosomes from the ATCC complete culture medium containing ED-FBS and purified them using the ExoQuick-TC PLUS kit. As shown in Figure 8A, exosomes are highly enriched in proteins in which tetraspanins (membrane proteins: CD63, CD9, CD81, CD82,

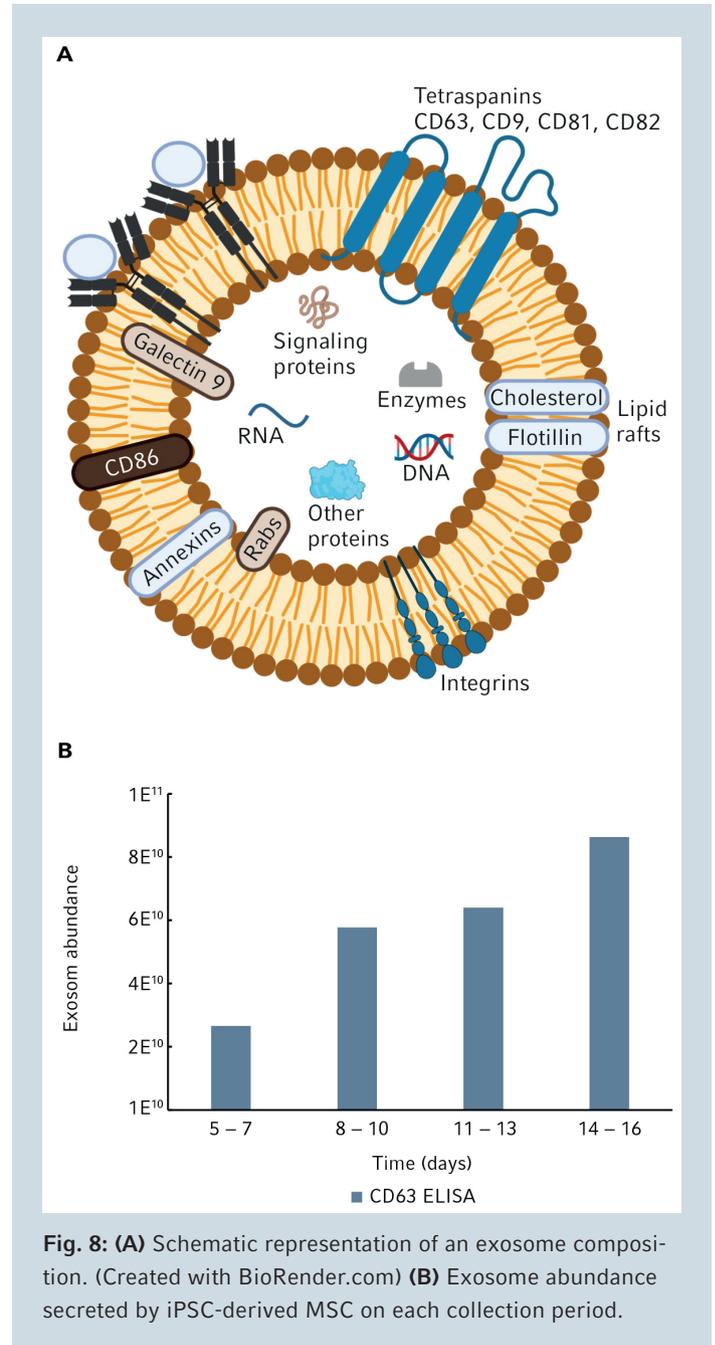


Fig. 8: (A) Schematic representation of an exosome composition. (Created with BioRender.com) (B) Exosome abundance secreted by iPSC-derived MSC on each collection period.

CD9, CD81, CD82) play a key role in cell invasion, penetration and fusion events. We used a direct Enzyme-Linked Immunosorbent Assay (ELISA), specifically the ExoELISA-ULTRA CD63 Kit to quantify the exosome abundance. We incubated the exosomes with a primary anti-CD63 antibody that binds the tetraspanin protein CD63 on the exosomal surface. The results show that the number of iPSC-derived MSC exosomes constantly increases from day 5 (2.6×10^{10}) to day 16 (8.6×10^{10} , Figure 8B). In addition, we found a direct correlation between cell density and secreted exosomes up to day 9 while the decrease in cell density in the bioreactor did not influence the exosome secretion after the stationary state.

To study the three-dimensional cell morphology on collagen coated microcarriers, we collected samples at days 5 and 9 of culture and visualized them through bright-field images (Figures 9A and B).

Taking into account the limitations of using bright-field systems for live cell imaging, we evaluated the morphology of human iPSC-derived MSC, demonstrating that at the early stage of the cell growth profile, a few cells attached on the microcarriers surface (Figure 9A) and progressively formed intermicrocarriers cellular bridges and aggregates (Figure 9B). It is well known that the 3D cellular environment allows the cells to mimic the *in vivo* cellular behavior, involving

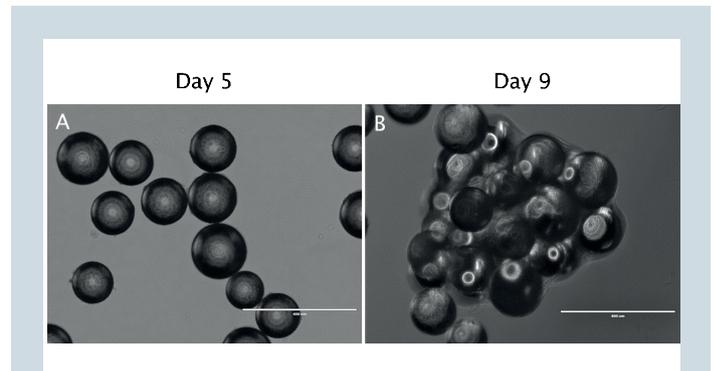


Fig. 9: Bright-field images (10x magnification) of human iPSC-derived MSC on collagen coated microcarriers at days 5 (A) and 9 (B) of culture.

cell-cell and cell-extracellular matrix (ECM) interactions, and promoting cell signaling and proliferation. It has recently been reported that although the mechanisms of exosome biogenesis, as well as the trafficking and release of vesicles are still not well understood, findings suggest that the exosome secretion and molecular cargoes can be altered by the cell microenvironment [12].

Conclusions

We have established the feasibility of producing MSC-derived exosomes using the SciVario twin bioreactor control system and BioBLU 1c Single-Use Bioreactors. With a working volume of up to 40 L and the possibility of further updates, the SciVario twin is a useful solution for your current and future bioprocessing needs. The design of the control system allows the precise manipulation of the cell culture environment, leading to rapid adhesion and proliferation of human

iPSC-derived MSC to the microcarrier's surface. We showed that by increasing the cell inoculation density, larger amounts of MSC exosomes were produced over time. These SciVario twin experiments are preliminary studies and have not yet been optimized to ascertain the maximum exosome production levels. However, our observations can serve as a guideline for further improvements in MSC-derived exosomes isolation, purification, and scale-up protocols.

Literature

- [1] Fernández O, Izquierdo G, Fernández V, et al. Adipose-derived mesenchymal stem cells (AdMSC) for the treatment of secondary-progressive multiple sclerosis: A triple blinded, placebo controlled, randomized phase I/II safety and feasibility study. *PLoS One* 2018; 13(5):e0195891.
- [2] Baboolal TG, et al. Synovial fluid hyaluronan mediates MSC attachment to cartilage, a potential novel mechanism contributing to cartilage repair in osteoarthritis using knee joint distraction. *Ann Rheum Dis* 2016; 75:908–915.
- [3] Wassef MA, et al. Therapeutic efficacy of differentiated versus undifferentiated mesenchymal stem cells in experimental type I diabetes in rat. *Biochemistry and Biophysics Reports* 2016; 5:468–475.
- [4] Cheng L, Zhao W, Hill AF. Exosomes and their role in the intercellular trafficking of normal and disease associated prion proteins. *Mol. Aspects Med.* 2018; 60:62–68.
- [5] Joo HS, Suh JH, Lee HJ, Bang ES, Lee JM. Current Knowledge and Future Perspectives on Mesenchymal Stem Cell-Derived Exosomes as a New Therapeutic Agent. *Int J Mol Sci.* 2020; 21:727.
- [6] Chen CC, Liu L, Ma F, et al. Elucidation of exosome migration across the blood–brain barrier model in vitro. *Cell Mol. Bioeng.* 2016; 9(4):509–529.
- [7] Yin K, Wang S, Zhao RC. Exosomes from mesenchymal stem/stromal cells: a new therapeutic paradigm. *Biomark Res.* 2019; 7:8.
- [8] June SH, Jinkwan K. Mesenchymal Stem Cell-derived Exosomes: Applications in Cell-free Therapy. *Korean J Clin. Lab. Sci.* 2018; 50:391-398.
- [9] Parallel Fed-batch CHO Culture on SciVario® twin, the Flexible Controller for All Your Bioprocess Needs. *Application Note* 2020; App-432.
- [10] Siddiquee K, Sha M. Microcarriers-based Expansion of Adipose-Derived Mesenchymal Stem Cells in Shake Flasks. *Bioprocessing Journal* 2013; 12:11-17.
- [11] Teixeira FG, Panchalingam KM, Anjo SI, et al. Do hypoxia/normoxia culturing conditions change the neuroregulatory profile of Wharton Jelly mesenchymal stem cell secretome? *Stem Cell Res Ther.* 2015; 6:133.
- [12] Thippabhotla S, Zhong C, He M. 3D cell culture stimulates the secretion of in vivo like extracellular vesicles. *Sci Rep* 2019; 9:13012.

Ordering information

Description	Order no. international
SciVario® twin Fermenter/Bioreactor Control System , base unit, 100 – 240 V/50/60 Hz, for 2 vessels	7600 100 001
BioBLU® 1c Single-Use Bioreactor , cell culture, open pipe, 2 pitched-blade impellers, optical pH, sterile, 4 pieces	1386 110 500
Adaptor , compression fitting, Pg 13.5 port to 12 mm sensor	M1287-5030
Centrifuge MiniSpin® plus , non-refrigerated, with Rotor F-45-12-11, 120 V/50 – 60 Hz (US)	0226 202 07*
Centrifuge 5430 R , keypad, refrigerated, with Rotor FA-45-30-11 incl. rotor lid, 120 V/50 – 60 Hz (US)	0226 206 01*
Incubator CellXpert® C170i , inner door with 4 door segments, handle right side, 1 – 20 % oxygen control, 100 – 127 V/50 – 60 Hz (JP/US)	6731 011 045*
Recirculation tube (from 3L decanter kit (M1287-119102))	M1287-9084
Adaptor , compression fitting, for Pg 13.5 port	M1287-5030
Eppendorf Conical Tubes , sterile, pyrogen-, DNase-, RNase-, human and bacterial DNA-free, colorless, 50 mL	0030 122 178
Eppendorf Conical Tubes , sterile, pyrogen-, DNase-, RNase-, human and bacterial DNA-free, colorless, 15 mL	0030 122 151
Eppendorf Cell Culture Flask, T-75 , TC treated, Filter cap, 80 pieces, 16 bags of 5 pieces	0030 711 122
Eppendorf Cell Culture Flask, T-175 , TC treated, Filter cap, 48 pieces, 12 bags of 4 pieces	0030 7121 29

*Inquire the part number for your country

Your local distributor: www.eppendorf.com/contact
 Eppendorf SE · Barkhausenweg 1 · 22339 · Hamburg · Germany
eppendorf@eppendorf.com · www.eppendorf.com/bioprocess

www.eppendorf.com/SciVario

