


Thermomixer comfort/Thermomixer R –Short Instructions

These short instructions may only be used by trained personnel who are familiar with the operating manual and the device. Observe the safety precautions.

Temperature Control and Mixing

Switch off "Progr." function (see below).

Use

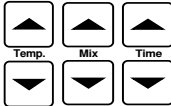

to set temperature, mixing frequency and/or time.




The device starts with temperature control/mixing. After the time has elapsed, mixing ends, however, temperature control continues.

Interval Mix

To switch on/off interval mix: Press  longer than 2 seconds.

Control light will light when function is activated.

Use

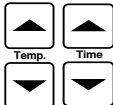

to set temperature and total time.



To set mix time



use the keys



To set stop time



The device starts alternates continually between temperature control and mixing. Mix – Pause, etc., until the total time has elapsed. Temperature control is then continued.


Thermomixer comfort/Thermomixer R – Short Instructions

Progr.

To switch on/off "Progr.": Press  longer than 2 seconds.

Control light will light when function is activated.

Use

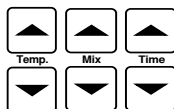


to set temperature, mixing frequency and/or time for the first mixing/temperature control run (temperature control and mixing or interval mix).



To switch to the second mixing/temperature control run.

Use

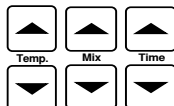

to set temperature, mixing frequency and/or time for the second mixing/temperature control run (routine run or interval mix).



To switch back to the first mixing/temperature control run.



The device starts with the first mixing/temperature control run. After the time has elapsed, it automatically starts with the second mixing/temperature control run. Temperature control is then continued with the temperature nominal value of the second run.

Short Mix



The device continues mixing with the set frequency as long as the key is pressed.

eppendorf

In touch with life

Your local distributor: www.eppendorf.com/worldwide

Eppendorf AG · Germany · Tel. +49 40 538 01-0 · E-Mail: eppendorf@eppendorf.com
Eppendorf North America, Inc. · USA · Toll free phone 800 645 3050 · E-Mail: info@eppendorf.com

Application Support

Europe, International: support@eppendorf.com · North America: support_NA@eppendorf.com
Asia, Pacific: support_AsiaPacific@eppendorf.com