Keep this rotor instruction manual and the centrifuge manual in the file.

Swing Rotors for Ultracentrifuge

P32ST2

Before using these rotors, please carefully read this instruction manual, "ROTORS, TUBES, BOTTLES AND CAPS", and the centrifuge instruction manual for its efficient operation and for your safety.

Keep this instruction manual for your reference and refer to it as required.

Contents

1. Specifications	_____	1
2. Applicable tubes	_____	2
3. How to use	_____	3
4. Maintenance	_____	6
5. Using cesium chloride	_____	7
6. Packing list	_____	8
7. Decontamination	_____	9
8. Rotor retirement	_____	9

Part No. [RLM rotor
9123190M

S998668-08/2025-10
Original instruction manual

Safety Reminder

- Centrifuge rotors rotating at high speed have considerable potential for damage to personal properties if used improperly.
For safe and proper use of this rotor, carefully read the centrifuge instruction manual and this rotor instruction manual before use and observe the instructions.

 **WARNING** : and  **CAUTION** : notes are used to call your attention in this manual to prevent personal injury or damage to the rotor and the centrifuge.
These notes are defined as follows.

-  **WARNING**: indicates a potentially hazardous situation which, if not avoided, could result in personal severe injury or possible death.
-  **CAUTION**: indicates a hazardous situation which, if not avoided, could result in personal injury or severe damage to the instrument.

WARNING

- Never use any material capable of producing flammable or explosive vapors. Your centrifuge is not explosion proof centrifuge and your rotor is not explosion proof rotor.
- Your centrifuge and rotor are not designed to confine any sample particles dispersed due to leakage. Therefore, when using toxic or radioactive samples or pathogenic or infectious blood samples, make sure to prepare necessary safety measures at your own responsibility.
- Never exceed the maximum speed of the rotor (mentioned on the rotor).
Always reduce rotor speed as instructed in this manual when rotor speed is limited due to sample density or kinds of tubes.
- Check the chemical resistance chart attached the rotor, and do not use any sample inapplicable to the rotor (including buckets). Using such a sample could corrode the rotor (including buckets).
- Do not remove the RLM adapter or optical adapter/disk from rotor, or replace it with the adapter/disk for another rotor. The adapter/disk is critical component that detects the over-speed of rotor. If an adapter/disk that is not compatible with the rotor is attached, the rotor could break, resulting in damage to the centrifuge.
- Do not allow the rotor temperature to rise 100°C.: This would cause the material brittle.
- Never sterilize the rotor by autoclaving or boiling. Otherwise the strength of rotor may be significantly decreased.
- When using swing rotor, be sure to set all buckets of the same type whether or not samples are put in buckets: Failure to do so could not only cause the rotor to vibrate, but could result in the rotor being deformed and the buckets becoming detached, which is very dangerous.
- If the centrifuge, rotor, or an accessory is contaminated by samples that are toxic or radioactive, or blood samples that are pathogenic or infectious, be sure to decontaminate the item according to good laboratory procedures and methods.
- If there is a possibility that the centrifuge, rotor, or an accessory is contaminated by samples that might impair human health (for example, samples that are toxic or radioactive, or blood samples that are pathogenic or infectious), it is your responsibility to sterilize or decontaminate the centrifuge, rotor, or the accessory properly before requesting repairs from an authorized sales or service representative.
- It is your responsibility to sterilize and/or decontaminate the centrifuge, rotor, or parts properly before returning them to an authorized sales or service representative.

CAUTION

- Check the chemical resistance chart attached to the rotor, and do not use any sample inapplicable to the tubes. Using such a sample could deteriorate them.
- The allowable speed may be lower depending on the density of sample, use of salt-such as cesium chloride-or the combination of tubes and tube adapters. Refer to this instruction manual.
- Always balance the tubes and samples within the allowable imbalance of rotor. Do not exceed the allowable imbalance. Otherwise the rotor and the centrifuge may be damaged.
- Do not load only one tube or load tubes asymmetrically: Asymmetrically loading may cause imbalance operation and damage the centrifuge and rotor.
- Clean the inside of the drive hole (crown hole) of the rotor and the surface of the drive shaft (crown) of the centrifuge once a month. If the drive hole or the drive shaft is stained or any foreign matter is adhered, the rotor may be improperly installed and come off during operation.
- Use the rotor tubes and bottles within their actual capacities.
- Do not use tubes that have exceeded their life expectancy. Failure to do so could result in damage of tubes and the rotor and the centrifuge. The life expectancy of tubes depends on factors such as the characteristics of samples, speed of the rotor used, and temperature. Always check for deterioration and damage (cracks, deformation, and so on) on tubes before using them. Do not use the tubes if you find such a problem.
- Inspect and maintain the rotor after use. If abnormality is observed, do not use it.
Contact an authorized sales or service representative.

1. Specifications

1.1 Construction

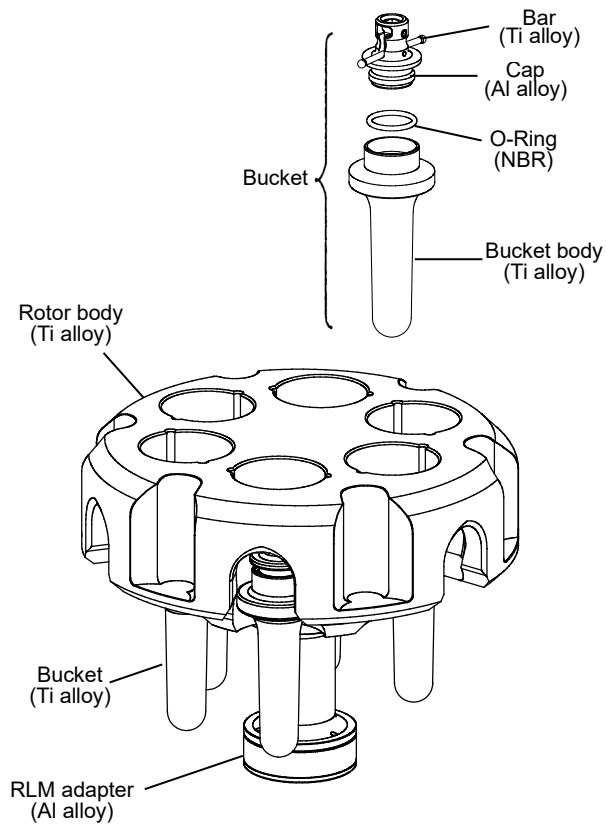
Fig. 1-1. Construction of P32ST2 rotor

1.2 Specifications

- Max. speed: 32,000 rpm
 - Max. RCF: 193,000 × g
 - Allowable density at maximum speed of the rotor: 1.2g/mL
 - K Factor at maximum speed of the rotor: 216
 - Max. rotor capacity(nominal): 16mL × 6tubes = 96mL
 - Material of the rotor body: Ti alloy
 - Material of the buckets: Ti alloy
 - Tube size: $\phi 16 \times 102$ mm
 - Weight(includes bucket and cap): 7.2 kg
 - (without bucket and cap) : 6.4 kg
 - Max. diameter: 21.7(34.0)cm
 - Max. height: 17.5cm
- (The numerical value in parentheses is maximum diameter measured with the bucket swinging.)

⚠ CAUTION:

- Do not use a RLM rotor with the centrifuge other than a series of a CP-NX/WX/MX ultracentrifuge.

1.3 Cross-sectional view of rotor

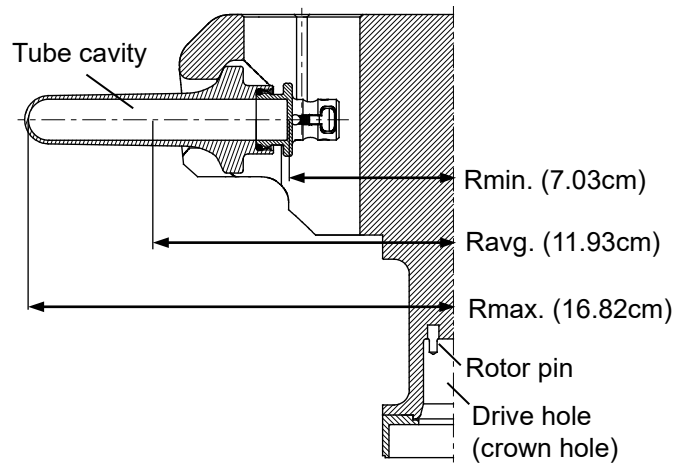
Fig. 1-2. Cross-sectional view of P32ST2 rotor

⚠ CAUTION:

Clean the inside of the drive hole (crown hole) of the rotor and the surface of the drive shaft (crown) of the centrifuge once a month. If the drive hole or the drive shaft is stained or any foreign matter is adhered, the rotor may be improperly installed and come off during operation.

1.4 Characteristic of rotor

Table 1-1. Characteristic of P32ST2

Rotor Speed (rpm)	RCF (×g)			K factor
	Rmax.	Ravg.	Rmin.	
5,000	4,700	3,330	1,960	8,843
10,000	18,800	13,300	7,860	2,211
20,000	75,200	53,300	31,400	553
24,000	108,000	76,800	45,300	384
30,000	169,000	120,000	70,700	246
32,000	193,000	137,000	80,500	216

※For detail about K factor, refer to "ROTORS, TUBES, BOTTLES AND CAPS(Part No.S999204)" attached to the centrifuge.

Calculation formula

$$RCF = 1.118 \times 10^{-5} \times R \times N^2 (\times g)$$

R: Rotating radius (cm)

N: Speed (rpm)

2. Applicable tubes

2.1 Applicable tubes

Table 2-1. Applicable tubes

Nominal capacity (mL)	Tubes				Tube adapter		Cap		Max. speed (rpm)	Max. RCF (xg)
	Part No.	Name	Size (φ X L cm)	Actual capacity (mL/tube)	Part No.	Name	Part No.	Name		
16	S303277A	16PET tube	1.6 X 10.2	17.0	---	---	---	---	32,000	193,000
16	332563A	16PA tube	1.6 X 10.2	15.0	---	---	---	---	32,000	193,000
15	S311533A	15PA seal tube	1.6 X 9.1	14.9	---	---	S413982	C6 space cap	32,000	193,000
3	336708A	3PA tube	1.1 X 4.7	2.8	336701A	3S16 adapter	---	---	24,500	77,500
4	329561A	4PA tube	1.1 X 6.0	3.7	336702A	4S16 adapter	---	---	24,800	86,800
5	332245A	5PA tube	1.3 X 5.2	4.9	336703A	5S16 adapter	---	---	24,800	80,900
3	S301599A	3PC tube	1.3 X 5.1	3.4	336703A	5S16 adapter	---	---	24,800	80,900
6.5	329445A	6.5PA tube	1.3 X 6.4	5.8	336704A	6.5S16 adapter	---	---	24,800	89,200
6.5	326415A	6.5PC tube	1.3 X 6.4	6.5	336704A	6.5S16 adapter	---	---	24,800	89,200
4.7	S303141A	4.7PC thick-walled tube	1.3 X 6.4	4.5	336704A	6.5S16 adapter	---	---	24,800	89,200

* Material of tube adapters: Polyacetal

* Material of the space cap: Aluminum alloy

1: When using 15PA seal tubes, we recommend the purchase of the tube kit, what is the set of tubes and necessary tools.

Refer to Chapter 6 for the tube kit.

CAUTION :


- Use the tubes and the adapters at a temperature between 4°C and 25°C. Otherwise, they may be broken or deformed during operation.
- Use the tubes that are same products, and the adapters that are same products. Otherwise, the imbalance operation may occur, and it may result in the damage to the centrifuge and the rotor.
- Chemical resistance and the strength of tubes varies with speed, temperature, and so on. Before using sample, fill tubes with water, buffer solution, etc. instead of sample and run them at the intended speed, temperature, and so on to ensure that there is no abnormality.

Refer to "ROTORS, TUBES, BOTTLES AND CAPS (Part No. S999204)" for the cleaning and sterilization of the tubes and the adapters.

How to use

How to use tubes/bottles varies in the kinds of tubes/bottles. The brief procedure of tubes/bottles handling is described as follows.

For the detail about handling tubes/bottles, refer to "ROTORS, TUBES, BOTTLES AND CAPS (Part No. S999204)" attached to the centrifuge.

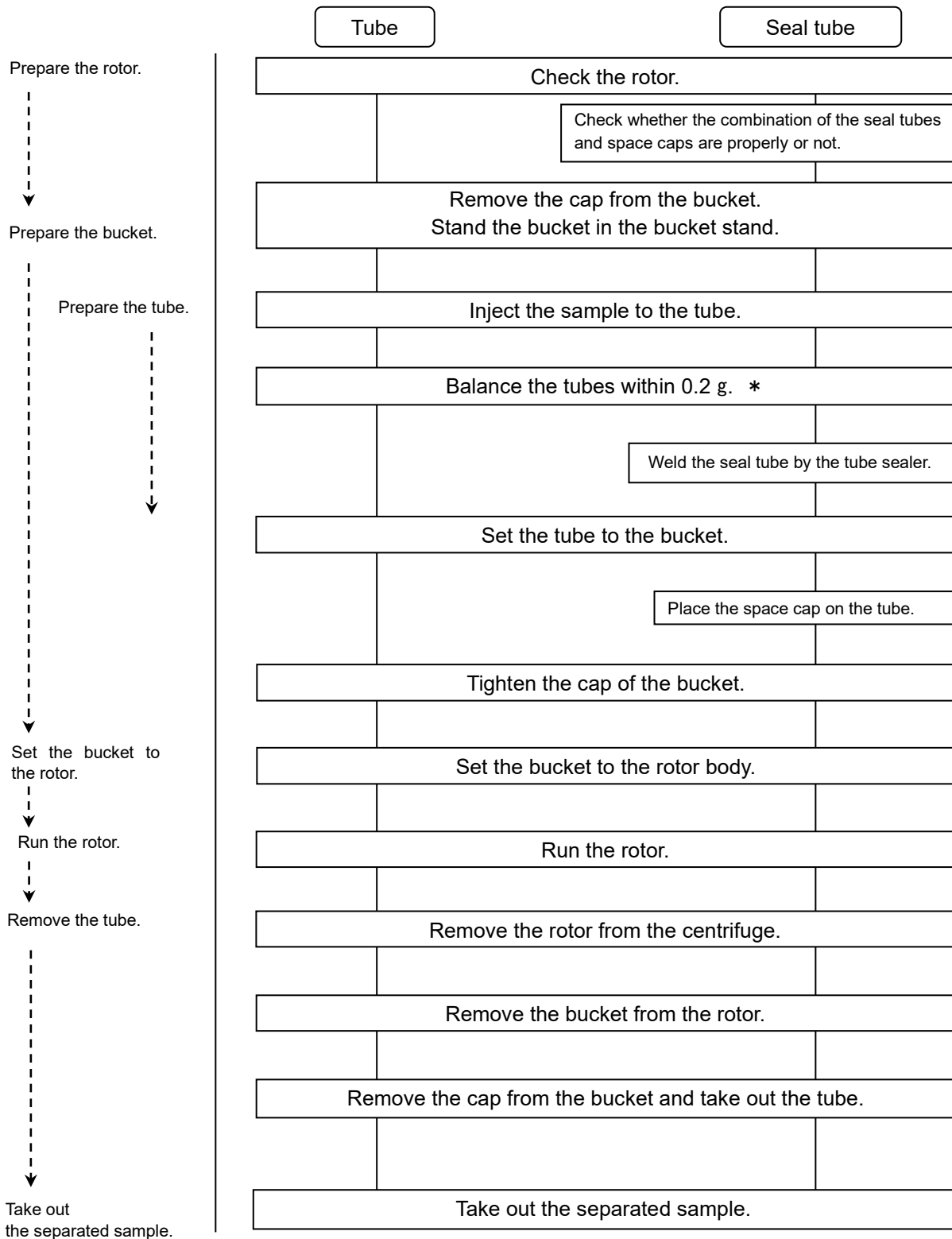(*) When using CP-NX series, CP-WX series, CP-MX series, approximately even quantities of sample in the tubes are sufficient for balancing, and extremely different sample quantities must be avoided. (The levels of samples in the tubes should be approximately equal. Their difference should be within 5mm.) For some kinds of tubes, if a small quantity of a sample is injected to them, the rotor can not be used at the allowable maximum speed. Tube deformation may occur due to a small quantity of a sample during operation.

=For handling=

● Allowable speed of rotors

1. To centrifuge a sample of an average density more than 1.2g/mL.
Determine the allowable speed from the following equation.

$$\text{Allowable speed (rpm)} = 32,000 \text{ (rpm)} \times \sqrt{\frac{1.2 \text{ (g/mL)}}{\text{Mean density of sample (g/mL)}}}$$

2. The allowable speed may be lower depending on the density of a sample when using of salt such as cesium chloride. Refer to "5. Using cesium chloride" for the detail.
3. The allowable speed may be lower depending on the combination of tubes and tube adapter, etc.
Refer to "2. Applicable tubes/bottles"
4. Reduce the rotor speed when the rotor has reached the primary life. Refer to "4. Maintenance ● Life of rotor".

● Tube/Bottle

Before operation, check whether the combination of tubes and tube adapters is properly or not.
Fill the tube up to within 5 mm from the top of the tube with sample.

● Seal tube

These tubes must be filled up the sample. Weld the seal tube completely.
Refer to "2. Applicable tubes" in this manual and "ROTORS, TUBES, BOTTLES AND CAPS (Part No. S99204)".

● Handling the rotor

WARNING :

- Never exceed the maximum speed of the rotor (mentioned on the rotor). Otherwise, rotor may be damaged and result in damage to the centrifuge.
- Check the chemical resistance chart attached the rotor (including buckets), and do not use any sample inapplicable to the rotor. Using such a sample could corrode the rotor (including buckets).

CAUTION :

- Inspect the rotor before use (see the "4. Maintenance · Inspecting rotor").
- Do not load only one tube or load tubes asymmetrically : Asymmetrically loading may cause imbalance operation and damage the centrifuge and rotor.
- Balance the tubes/bottles and the samples including the caps and the tube adapters within the allowable imbalance of the rotor.
Do not exceed the allowable imbalance. Otherwise, the centrifuge or the rotor may be damaged due to excessive imbalance operation.
- If foreign matter adheres to the rotor (including buckets), wipe it off promptly. Otherwise, the rotor (including buckets) may be corroded.


Set the balanced tubes to the rotor symmetrically.

<Tightening the cap to the bucket>

CAUTION :

Imperfect tightening of the cap may cause imbalance operation and it may happen the serious accident.

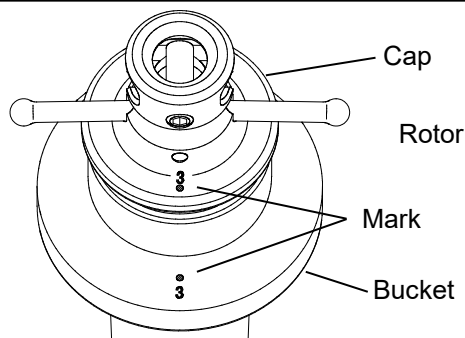Tighten the numbered cap to the same numbered bucket to match each mark of the bucket and the cap.


Tighten the cap until the bottom of the cap bumps to the bucket and the cap must be adjusted to the match mark of the bucket.

<Setting the bucket to the rotor>

CAUTION :

Check that "P32ST2" is marked on the buckets, caps, and the rotor before use.

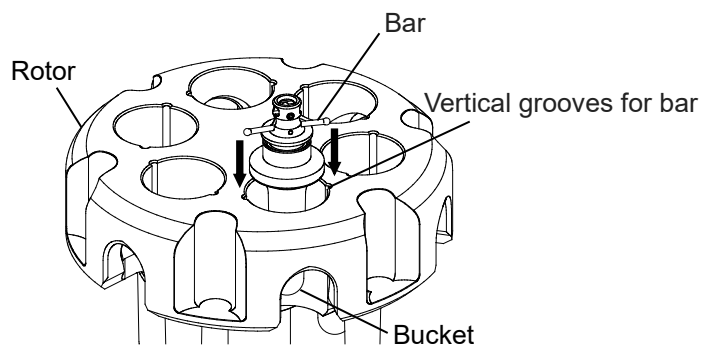
Insert a bar of the bucket to vertical groove of the rotor.

Match the numbered buckets with number of rotor and set the buckets to the rotor.

⚠ WARNING :

- Be sure to set all buckets of the same type whether or not samples are put in buckets: Failure to do so could not only cause the rotor to vibrate, but could result in the rotor being deformed and the buckets becoming detached, which is very dangerous. Never use buckets made by other companies, or any other type of bucket that is not exclusively made for the rotor even if it is made by us.
- Make sure to use the buckets with the cap marked "P32ST2" which are designed exclusively for the P32ST2 rotors. Otherwise, the rotor and buckets may be damaged.

<Examples of incorrect bucket setting>

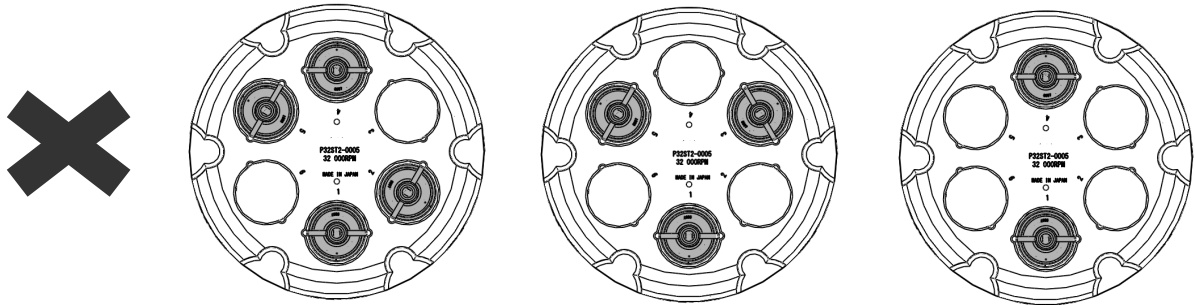Tighten the cap to the bucket securely. Set the bucket to the rotor securely.

⚠ WARNING :

Before setting a rotor in ultracentrifuge, make sure that buckets have been securely installed in the rotor: Incorrect installation of buckets in rotor could damage the rotor or detach the buckets, which is very dangerous.

⚠ CAUTION :

Inspect and maintain the rotor after use. If abnormality is observed, do not use it. Contact an authorized sales or service representative.

● **Sample limitation**

⚠ WARNING :

- Never use any material capable of producing flammable or explosive vapors. Your centrifuge and rotor are not designed to confine any sample particles dispersed due to leakage.
- When using toxic or radioactive samples or pathogenic or infectious blood samples, make sure to prepare necessary safety measures at your own responsibility.

● **RLM adapter**

⚠ WARNING :

Do not remove the RLM adapter from rotor, or replace it with the adapter/disk for another rotor. The adapter is a critical component that detects the over-speed of rotor: If an adapter that is compatible with the rotor is attached, the rotor could break, resulting in damage to the ultracentrifuge.

⚠ CAUTION :


- Do not bring the RLM adapter near a magnet, or scratch it: Doing so will erase the memory stored in adapter, and make the rotor unusable. To prevent the adapter from being scratched, store the rotor with RLM adapter, using the stand provided with the rotor (rotor stand for protecting adapter). If the adapter is scratched, immediately stop using rotor and contact an authorized sales or service representative.
- If sample etc. adheres to the RLM adapter of rotor, immediately wipe it off with a soft cloth etc. If the RLM adapter corroded, the rotor will be unusable.

4. Maintenance

● Maintenance of rotor

⚠ WARNING :
Do not allow the rotor (including buckets) temperature to rise over 100°C. Otherwise, the rotor can be brittle.

⚠ CAUTION :
Use a neutral detergent having a pH between 5 and 9. Otherwise, the rotor can be corroded.

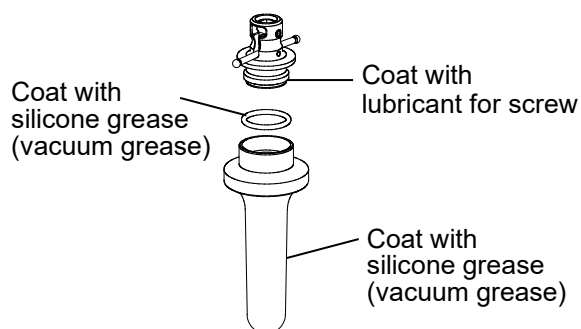Wash the bucket with tap water or a dilute solution of neutral detergent and rinse it out with distilled water. Turn the bucket upside down with the cap removed to dry it well. Check the bucket is completely dried and then put the light coat of silicone grease on the surface of the bucket, the light coat of lubricant for screw on the thread portion, the light coat of silicone grease (vacuum grease) on the O-ring and the gasket. Store them in a dry place. If foreign matter adheres to the bucket, soak the bucket in warm water in one or two hours and remove foreign matter with a soft brush etc. Use a neutral detergent having a pH between 5 and 9. If you can not remove it, contact an authorized sales or service representative.

Wipe the stain off the rotor body and then put the light coat of silicone grease (vacuum grease) on the surface of the rotor body. Store it in a dry place.

Put the light coat of silicone grease inside the drive hole (crown hole) and then wipe it out with a soft cloth.

● Sterilizing rotor

⚠ WARNING :
Never sterilize the rotor (including buckets) by autoclaving or boiling. Otherwise the strength of rotor (including buckets) may be significantly decreased.

Sterilize the rotor according to either gas sterilization method (ethylene oxide or formaldehyde), chemical sterilization method (70% ethanol, 3% hydrogen peroxide, 3% formalin), or ultraviolet rays sterilization.

* Do not dip the rotor in the formalin (3%) solution more than 2 hours.

● Inspecting rotor

Check the rotor after every run as follows.

- The thread portion is worn. Contact an authorized sales or service representative.
- The tube cavity is corroded. Contact an authorized sales or service representative.
- The O-ring, the gasket, etc. are deteriorated. Replace it with the new one.
- The adapter is stained. Wipe it out with a soft cloth.

Check the rotor every 100 hours as follows.

- The rotor pin is bent. Contact an authorized sales or service representative.

Check the rotor once a month as follows..

- The inside of drive hole (crown hole) is stained. Wipe it out with a soft cloth dampened with water.

⚠ CAUTION :
Clean the inside of the drive hole (crown hole) of the rotor and the surface of the drive shaft (crown) of the centrifuge once a month. If the drive hole or the drive shaft is stained or any foreign matter is adhered, the rotor may be improperly installed and come off during operation.

● Life of rotor

While using rotor repeatedly, it's strength decreases gradually due to fatigue and creep of material, by which the rotor life is determined.

Sum up the number of operations and number of operating hours. And if either value reaches the figures (the primary life) shown in Table 4, decrease the maximum speed to 28,000 rpm. When the rotor has reached the primary life, maintenance of the rotor is required at extra cost (extra-cost option).

If the operations or operating hours of rotor reach the figures in Table 4 after decreasing the maximum speed to 28,000 rpm (Secondary life), don't use the rotor any more. If you use the RLM rotor at a speed lower than its maximum permissible speed, the life of the RLM rotor will be prolonged automatically by calculating from the strength of the rotor material.

Table 4. Rotor life

500 operations	2,500 hr
----------------	----------

5.Using cesium chloride

The cesium chloride (CsCl) solution is frequently used as a density gradient medium, but if the CsCl solution with a high density is used, it may saturate during the rotation of the rotor depending on the rotation condition which may cause CsCl crystals to be educed. Educated crystals have a high density(approx.4g/mL) and apply an excessive load to the rotor which is very dangerous.

The crystallization will also greatly affect the density gradient of the CsCl solution and separation status of the sample. Therefore, always use a solution which will not crystallize the allowable speed.

Fig .5-1 shows the relationship between the maximum allowable density and the speed at which cesium chloride does not crystallize in the condition that standard tubes is filled with the cesium chloride solution. In this diagram, crystallization will occur above the curves.

Be sure to select the density of the cesium chloride solution and the rotor speed at which cesium chloride does not crystallize. The curve in this diagram vary in temperature etc. Before actual operation, centrifuge the cesium chloride solution with the same density to ensure that crystallization will not occur.

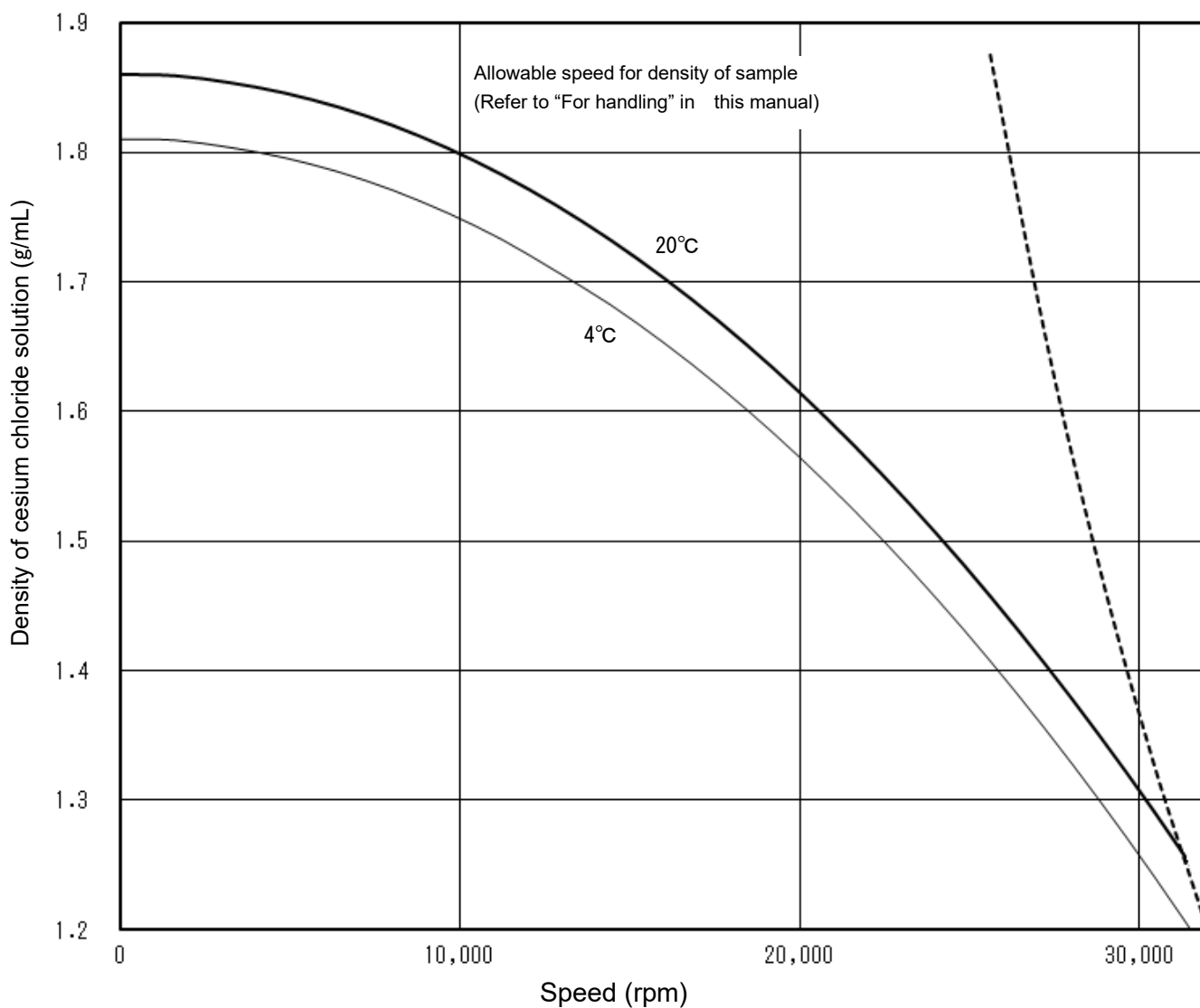
Fig.5-1. Relationship between density of cesium chloride and speed

6. Packing list

Check the parts and accessories of rotor when receiving the rotor.
 Contact an authorized sales or service representative when the need arises.

No.	Name	Qty	Remark	Part No. (Qty)	No.	Name	Qty	Remark	Part No. (Qty)
①	P32ST2 Rotor	1			⑧	Rotor stand	1		216715 (1)
②	Rotor body	1			⑨	Tool box(C)	1		
③	Bucket body	6*			⑩	Tweezers	1		452701 (1)
④	NBR O-Ring	6*		8062020 (1)	⑪	Bucket stand	1		209642A (1)
⑤	Cap	6*			⑫	Lubricant for screw	1		84810601 (1)
⑥	NBR O-Ring	18	For spare		⑬	Silicone grease (vacuum grease)	1		483719 (1)
⑦	16PA tube	2	50pcs./set	332563A (50)	⑭	Instruction manual	1		S998668 (1)


*: the quantities of the parts that are mounted to the rotor

When using 15PA Seal Tubes, we recommend the purchase of under-mentioned tube kit, what is the set of tubes and necessary tools.

● 15PA Seal Tube Kit (Part No.S311535A)
 (component)

S311533A 15PA Seal Tube : 2 (50 pcs./pack)
 S413982 C6 Space Cap : 6
 S407157 Tube Setter : 1

※When using 15PA Seal Tubes, Tube Sealer (STF3, STF2) and Tube Rack (C3) (Part No.S206017F) are needed.

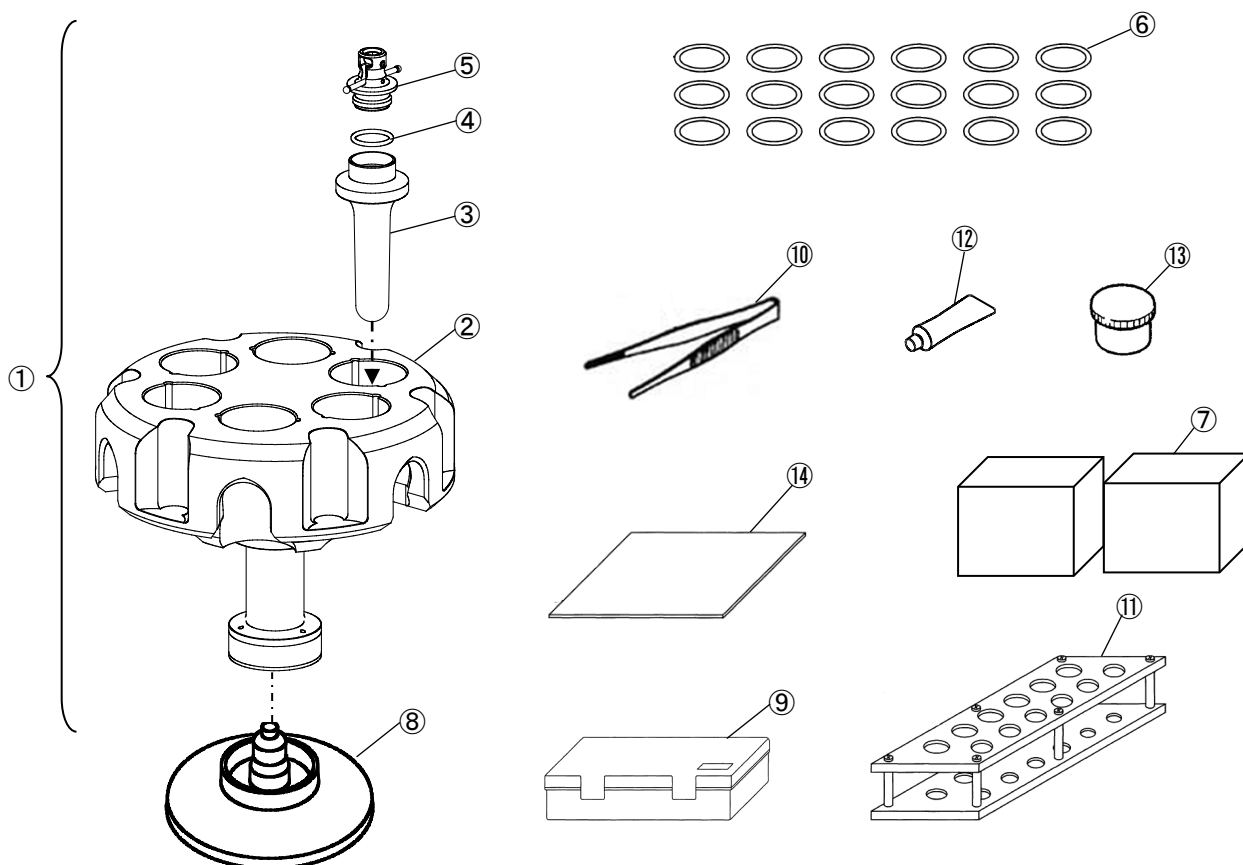7. Decontamination

⚠ WARNING:

- If the centrifuge, rotor or an accessory is contaminated by samples that are toxic or radioactive, or blood samples that are pathogenic or infectious, be sure to decontaminate the item according to good laboratory procedures and methods.
- If there is a possibility that the rotor or an accessory is contaminated by samples that might impair human health (for example, samples that are toxic or radioactive, or blood samples that are pathogenic or infectious), it is your responsibility to sterilize or decontaminate the rotor or the accessory properly before requesting repairs from an authorized sales or service representative. Note that we cannot repair the centrifuge, rotor or the accessory unless sterilization or decontamination is completed.
- It is your responsibility to sterilize and/or decontaminate the rotor or parts properly before returning them to an authorized sales or service representative. In such cases, copy the attached decontamination sheet and fill out the copied sheet, then attach it to the item to be returned. We may ask you about the treatment for the rotor or parts if the decontamination is checked and judged as insufficient by us. It is your responsibility to bear the cost of sterilization or decontamination. Note that we cannot repair or inspect the rotor or the accessory unless sterilization or decontamination is completed.

8. Rotor retirement

After many years of use, there will be inevitably some corrosion or stress corrosion. At some points, the combination of such damage and metal fatigue could make the rotor vulnerable to a failure. Although a rotor may appear to be in a good condition, you should follow the rotor retirement recommendation shown below.

Rotor	Material	Retire After Years
Angle rotor	Titanium alloy	12
Vertical rotor	Aluminum alloy	10
Swing rotor	Titanium alloy Aluminum alloy	10
Zonal rotor	Titanium alloy	10
Continuous flow rotor	Titanium alloy	10

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Copyright © 2025 Eppendorf Himac Technologies Co., Ltd. All rights reserved.

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