

Register your instrument!
www.eppendorf.com/myeppendorf

Rotor FA-10x50 ID

Instructions for Use

Copyright © 2024 Eppendorf SE, Germany. All rights reserved, including graphics and images. No part of this publication may be reproduced without the prior permission of the copyright owner.

Eppendorf® and the Eppendorf Brand Design are registered trademarks of Eppendorf SE, Germany.

Eppendorf trademarks and trademarks of third parties may appear in this manual. All trademarks are the property of their respective owners. The respective trademark name, representations and listed owners can be found on www.eppendorf.com/ip.

U.S. Patents and U.S. Design Patents are listed on www.eppendorf.com/ip.

Table of contents

1	About this manual	5
1.1	Notes on this manual.....	5
1.2	Warning notice structure.....	5
1.3	Graphics.....	5
1.4	Other applicable documents.....	6
2	Safety	6
2.1	Personal injury.....	6
2.1.1	Biological hazards.....	6
2.1.2	Mechanical hazards.....	6
2.2	Material damage.....	6
2.2.1	Chemical hazards.....	6
2.2.2	Mechanical hazards.....	7
2.2.3	Incorrect handling.....	7
2.2.4	Loss of samples.....	7
2.3	Information on product liability.....	7
3	Delivery package	8
4	Product description	8
4.1	Features.....	8
4.2	Product overview.....	8
5	Application	9
5.1	k-factor.....	9
5.2	Rotor table.....	9
5.3	Selecting the rotor and accessories.....	12
5.4	Inserting the rotor.....	12
5.5	Removing the rotor.....	13
5.6	Opening the rotor lid.....	13
5.7	Closing the rotor lid.....	13
5.8	Loading the rotor.....	13
6	Maintenance	14
6.1	Maintenance plan.....	14
6.2	Maintenance.....	14
6.2.1	Checking the rotor and accessories.....	14
6.3	Cleaning.....	15
6.3.1	Cleaning a contaminated fixed-angle rotor.....	15
6.3.2	Cleaning a fixed-angle rotor after tube breakage.....	15
6.4	Decontamination.....	16
6.4.1	Decontaminating the rotor.....	16
6.4.2	Autoclaving the rotor.....	17
7	Transport	17
7.1	Shipping the rotor.....	17

8 Disposal	18
8.1 Legal regulations.....	18
8.2 Preparing for disposal.....	19
9 Technical data	19
9.1 Application parameters.....	19
9.2 Service life.....	19
10 Ordering information	20

1 About this manual

1.1 Notes on this manual

This document describes the rotor and does not replace the operating manual for the device.

1. Before using the rotor, read these instructions for use.
2. Also read the operating manual for the device with which you are using the rotor.



Eppendorf and third-party trademarks may appear in these instructions for use. All trademarks are the property of their respective owners. The respective trademarks, figurative marks and listed owners can be found at www.eppendorf.com/ip.

1.2 Warning notice structure



HAZARD LEVEL! Type of danger

- Source of danger
- Consequences of disregarding the danger
- Measures to avoid the danger

Symbol	Hazard level	Type of danger	Meaning
	DANGER	Personal injury	Will lead to severe injuries or death.
	WARNING	Personal injury	May lead to severe injury or death.
	CAUTION	Personal injury	May lead to minor or moderate injuries.
	NOTE	Material damage	May lead to material damage.

1.3 Graphics

Depiction	Meaning
1.	Work steps
2.	
•	Bullet point
<i>Text</i>	Display text
Key	Name for port, button, status lamp, or key
	Important information
	Hint

1.4 Other applicable documents

The following documents supplement these instructions for use:

- Operating manual for SpinPro 6 R
- Operating manual for SpinPro 7 R

2 Safety

2.1 Personal injury

2.1.1 Biological hazards

If unsuitable or damaged rotors, rotor lids or accessories are used in aerosol-tight centrifugation, hazardous aerosols may escape.

- Only use aerosol-tight rotors and rotor lids and aerosol-tight accessories.
- Check all rotors, rotor lids and accessories for perfect condition before each use.

2.1.2 Mechanical hazards

If you handle rotors incorrectly, the rotors may fall and injure people.

- Hold rotors with both hands.

If rotors break, sharp-edged fragments can be created that can injure people.

- Wear cut-resistant gloves when collecting sharp-edged fragments.
- Avoid material damage to the rotors.

2.2 Material damage

2.2.1 Chemical hazards

Aggressive chemicals as well as aggressive cleaning agents and disinfectants can damage the material of the device, rotor and accessories. This may cause damage to the interior of the device during operation.

- Protect all components from aggressive chemicals.
- Clean and disinfect all components only with recommended cleaning agents and disinfectants.
- Clean the rotor, device and accessories immediately if they come into contact with aggressive chemicals.
- Check all components for material changes before each use.
- Take a damaged device out of operation.
- Replace a damaged rotor.
- Replace damaged accessories.

2.2.2 Mechanical hazards

A damaged rotor or damaged accessories can be destroyed when used in the device.

- Check the rotor and accessories for visible damage and material changes before each use.
- Use only undamaged rotors and accessories.
- Use rotors and accessories only within their service life.

2.2.3 Incorrect handling

Improperly attached rotors can damage the device, rotor and tubes.

- Insert the rotor as described in this manual.
- Tighten the rotor nut.
- Close the rotor lids and caps, if present, as described in this manual.
- Stop centrifugation immediately if you hear any unusual noises.

2.2.4 Loss of samples

Incorrect handling of tubes may result in broken glass and sample loss.

- Use suitable tubes.
- Insert the tubes with an adapter as necessary.

Incorrect handling of tubes may cause sample loss.

- Observe the manufacturer's instructions for tubes.
- For centrifugation, only use liquids with a maximum density of 1.2 g/mL.

2.3 Information on product liability

The owner of the device will be held liable for personal and material damage in the following cases:

- The device is used outside of its intended use
- The device is not used in accordance with the operating manual
- Manipulation of safety devices
- The device has spare parts installed that are not authorized by Eppendorf SE
- The device is used with accessories or consumables that are not recommended by Eppendorf SE
- Cleaning agents are used that are not recommended by Eppendorf SE
- Chemicals are used that are not recommended by Eppendorf SE
- Shipment not in original packing or in improper substitute packing
- The device is maintained or repaired by persons not authorized by Eppendorf SE
- Unauthorized modifications

3 Delivery package

Quantity	Description
1	Rotor with lid
1	Pivot grease – 5 mL
1	Instructions for use

4 Product description

4.1 Features

The FA-10×50 ID rotor is a fixed-angle rotor with a capacity of up to 50 mL per bore.

The rotor has the following specific features:

- High-speed rotor
- Rotor identification by RFID
- Aerosol-tight centrifugation
- Eppendorf QuickLock Pro rotor lid lock

The following consumables can be used with the rotor:

- 9 mL tubes
- 5 mL – 50 mL conical tubes
- 2.6 mL – 30 mL round-bottom tubes

4.2 Product overview

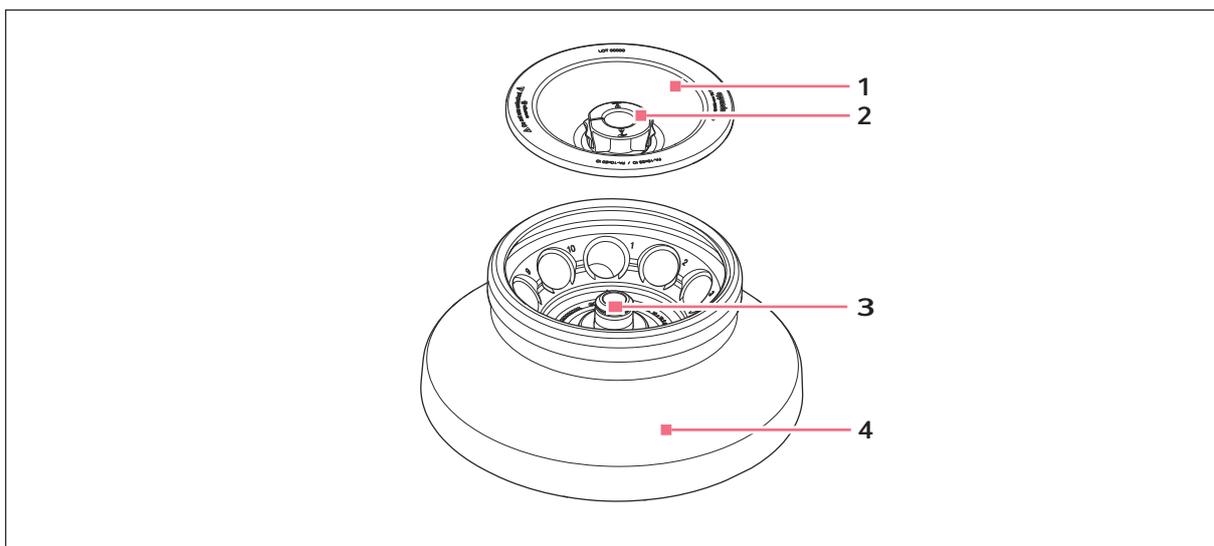


Fig. 4-1: Components of rotor FA-10×50 ID

1 Rotor lid

3 Rotor hub

2 Unlocking mechanism

4 Rotor body

5 Application

5.1 k-factor

The k-factor is a rotor-specific factor for determining the sedimentation efficiency and depends on the position of the reaction tube in the rotor. The k-factor takes into account the g force and the distance. The smaller the k-factor, the shorter the centrifugation time.

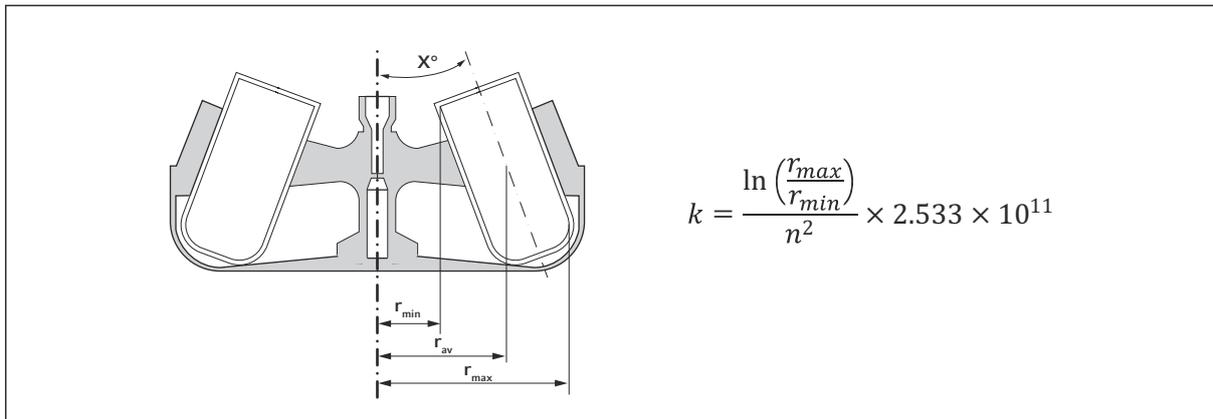


Fig. 5-1: Determination of the k-factor

r_{max} Maximum centrifugation radius

α Rotation angle

r_{min} Minimum centrifugation radius

n Maximum rotation speed in rpm


r_{av} Average centrifugation radius





The k-factors can be found in the rotor tables.

5.2 Rotor table

Suitable tubes for the FA-10x50 ID rotor

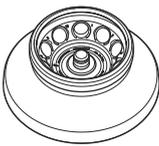FA-10x50 ID		
	Max. g -force	$25012 \times g$
	Max. speed	12420 rpm
	Max. load (adapter, tube and contents)	10×74.5 g

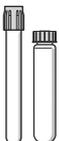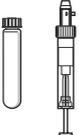Tube	Tube Capacity Tubes per adapter/rotor	Adapter Order no. (International)	Bottom shape Tube diameter Max. tube length with tube lid	Max. <i>g</i> -force Radius	k-factor
	Round-bottom tube 16 mL 1/10	 5820 720 000	Round Ø 18.1 mm 120 mm	24253 × <i>g</i> 14.1 cm	1157
	Round-bottom tube Ø 13 mm × 75 mm) 2.6 mL – 5 mL 1/10	 5820 726 008	Round Ø 13.5 mm 125 mm	24322 × <i>g</i> 14.1 cm	292
	Round-bottom tube Ø 13 mm × 100 mm 4 mL – 8 mL 1/10	 5820 725 001	Round Ø 13.5 mm 125 mm	24322 × <i>g</i> 14.1 cm	438
	Reaction vessel 5 mL 1/10	 5820 730 005	Conical Ø 17 mm 125 mm	24322 × <i>g</i> 14.1 cm	450
	Round-bottom tube Ø 16 mm × 75 mm 5.5 mL – 10 mL 1/10	 5820 728 000	Round Ø 16 mm —	24322 × <i>g</i> 14.1 cm	393
	Round-bottom tube Ø 16 mm × 100 mm 7.5 mL – 12 mL 1/10	 5820 727 004	Round Ø 16 mm 125 mm	24322 × <i>g</i> 14.1 cm	532

Tube	Tube Capacity Tubes per adapter/rotor	Adapter Order no. (International)	Bottom shape Tube diameter Max. tube length with tube lid	Max. <i>g</i> -force Radius	k-factor
	Tube 9 mL 1/10	 5820 729 007	Round Ø 17.5 mm 125 mm	24322 × <i>g</i> 14.1 cm	438
	Conical tube 15 mL 1/10	 5820 717 009	Conical Ø 17 mm 130 mm	24650 × <i>g</i> 14.3 cm	1332
	Round-bottom tube 30 mL 1/10	 5820 721 006	Round Ø 25.7 mm 104 mm	21769 × <i>g</i> 12.6 cm	963
	Conical tube 35 mL 1/10	 5820 722 002	Conical Ø 28.7 mm 113 mm	22770 × <i>g</i> 13.2 cm	1060
	Conical tube 25 mL 1/10	 5820 734 000	Conical Ø 29.8 mm 88 mm	20458 × <i>g</i> 11.9 cm	1039
	Conical tube 25 mL 1/10	 5820 733 004	Conical Ø 29.8 mm 83 mm	21683 × <i>g</i> 12.6 cm	676
	Conical tube 50 mL 1/10	—	Conical Ø 29.6 mm 125 mm	25012 × <i>g</i> 14.5 cm	1138
	Conical tube 50 mL 1/10	—	Conical Ø 30 mm 125 mm	25012 × <i>g</i> 14.5 cm	1083

5.3 Selecting the rotor and accessories

1. Select the rotor, adapter and tubes according to your application.
2. Use a rotor lid suitable for the rotor. Ensure that the rotor designation on the rotor and the rotor lid match.
3. Ensure that rotors, the rotor lid and accessories meet the following requirements:
 - Maximum service life not exceeded
 - Intact
 - No corrosion
 - No material changes
 - Undamaged and clean seals
 - Seals not protruding at any point
4. Make sure the tubes comply with the following requirements:
 - Approved by the manufacturer for this application
 - Chemically resistant to the sample
 - Designed to withstand the intended load
 - Intact
 - No deformation
 - Not brittle

5.4 Inserting the rotor

The rotor is detected automatically when inserted, it is equipped with RFID rotor identification. The device limits the speed of the detected rotor to its permissible maximum speed. The rotor detection allows the rotor to count the number of cycles used. The rotor detection prevents a rotor from being used in a device that is incompatible with the rotor. The rotor also prevents use after the number of cycles has been exceeded.

Tool:

- Rotor key

Prerequisites:

- A usable rotor is available.

1. Place the rotor vertically on the motor shaft.

The rotor is detected.

2. Tighten the rotor nut. To do this, turn the rotor key clockwise.

5.5 Removing the rotor


Tool:

- Rotor key

1. Loosen the rotor nut by turning the rotor key counterclockwise.
2. Lift the rotor vertically out of the device.

5.6 Opening the rotor lid


Opening the rotor lid with Eppendorf QuickLock Pro

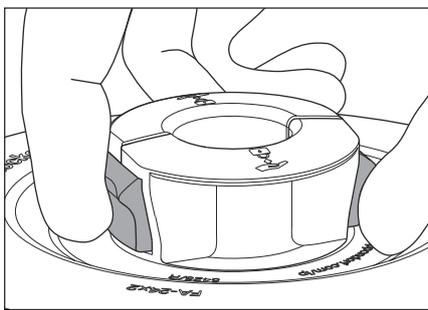

1. Press the locking mechanism on both sides simultaneously.
2. Remove the rotor lid.

5.7 Closing the rotor lid

Closing the rotor lid with Eppendorf QuickLock Pro

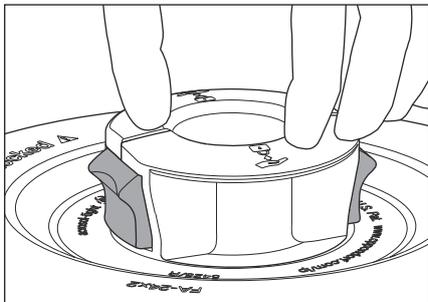

1. Ensure that the rotor designation on the rotor and the rotor lid match.
2. Place the rotor lid on the rotor.
3. Push the handle down until the rotor lid audibly engages and the unlocking mechanism is back in its original position.

5.8 Loading the rotor



When using different tube types, the weight the rotor is loaded with must be distributed symmetrically.



NOTICE! Sample loss

Open tube lids can break off during centrifugation.

- Before centrifugation, make sure that the tube lids are closed correctly.


**NOTICE! Damage to device**

If rotor are loaded incorrectly, the rotor and the device may become damaged during centrifugation.

- Load the rotor symmetrically.
- Only load the rotor to its maximum load capacity.

Prerequisites:

- The tubes are filled.
- The filled tubes and adapters have the same weight.

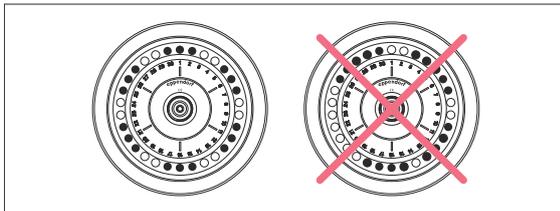

1. Place the tubes in suitable adapters if necessary.
2. Insert tubes of the same weight in the rotor holes correctly aligned and symmetrically.

6 Maintenance

6.1 Maintenance plan

Interval	Maintenance work
As required	☞ Chapter 6.3.1 "Cleaning a contaminated fixed-angle rotor" on page 15
	☞ Chapter 6.3.2 "Cleaning a fixed-angle rotor after tube breakage" on page 15
	☞ Chapter 6.4.1 "Decontaminating the rotor" on page 16
	☞ Chapter 6.4.2 "Autoclaving the rotor" on page 17
Before each use	☞ Chapter 6.2.1 "Checking the rotor and accessories" on page 14

6.2 Maintenance

6.2.1 Checking the rotor and accessories

1. Check the rotors and accessories for damage.
 - Take damaged rotors and damaged accessories out of operation.
2. Check the service life of the rotors and the accessories.
 - Replace rotors and accessories that have exceeded their service life.

6.3 Cleaning

6.3.1 Cleaning a contaminated fixed-angle rotor



Clean the rotor by hand.

Do not submerge the rotor under water. When submerged, liquid can enter the gaps of the rotor.

Material:

- Deionized water
- Soap-based cleaning agent
- Cloth
- Cleaning brush

Prerequisites:

- The rotor has been removed from the device.

1. Wet the cloth and cleaning brush with deionized water and cleaning agent.
2. Clean the rotor body and rotor lid.
3. To remove the cleaning agent from the rotor body and rotor lid, wipe these off with deionized water.
4. Let the rotor body and rotor lid dry.



Cleaning can dry out the seal areas on the rotor. This can result in sluggish rotor lid operation. If this is the case, lightly grease the seal with the supplied pivot grease.

6.3.2 Cleaning a fixed-angle rotor after tube breakage



Clean the rotor by hand.

Do not submerge the rotor under water. When submerged, liquid can enter the gaps of the rotor.



Very small glass particles can deposit in adapters and parts made of rubber, e.g., seals. Clean these parts very carefully.



WARNING! Contamination

Biological substances leak from broken tubes. People who come into contact with these substances may suffer damage to their health.

- Wear your personal protective equipment.
- If necessary, carry out decontamination before cleaning.

Material:

- Water
- Soap-based cleaning agent
- Cloth

1. Remove all glass splinters.
2. Remove the rotor and the accessories.
3. To remove any shards from tubes, clean the rotor, the rotor bores and the accessories with a damp cloth.
4. Replace the adapters to prevent any further damage.
5. Allow the rotor and the accessories to dry.



Cleaning can dry out the seal areas on the rotor. This can result in sluggish rotor lid operation. If this is the case, lightly grease the seal with the supplied pivot grease.

6.4 Decontamination

6.4.1 Decontaminating the rotor



Decontamination with at least 70% ethanol is described here. Observe the manufacturer's information when using an alternative disinfectant.

Material:

- Disinfectant with at least 70 % ethanol
- Deionized water
- Cloth

Prerequisites:

- The rotor has been removed from the device.

1. Wet the cloth with disinfectant.
2. Wipe the inside of the rotor and the accessories with the cloth.
3. Wipe off the disinfectant with deionized water.
4. Clean the rotor and the accessories.



Cleaning can dry out the seal areas on the rotor. This can result in sluggish rotor lid operation. If this is the case, lightly grease the seal with the supplied pivot grease.

6.4.2 Autoclaving the rotor

The following are autoclavable:

- Rotor body
- Rotor lid with seal

Protective equipment:

- Laboratory protective clothing

1. Autoclave the rotor body and the rotor lid at 121 °C with steam for 20 min.



- You can also autoclave the rotor body at 134 °C with steam for 20 min.
- Cleaning can dry out the seal areas on the rotor. This can result in sluggish rotor lid operation. If this is the case, lightly grease the seal with the supplied pivot grease.

7 Transport

7.1 Shipping the rotor



Use the original packing for transport. In the event that the original packing is no longer available, ensure that the rotor is adequately protected by replacement packing during storage and further transport. Eppendorf SE is not liable for damage caused by improper replacement packing.



WARNING! Contamination

If you store or ship a contaminated rotor, people may be contaminated and suffer damage to their health.

- Decontaminate the rotor before storage or shipping.

Material:

- Packing

Prerequisites:

- The rotor has been decontaminated and cleaned.

1. Download the decontamination certificate for return shipments from the webpage www.eppendorf.com.
2. Complete the decontamination certificate.
3. Pack the rotor.
4. Attach the decontamination certificate to the outside of the packing in a way that is safe for transport.
5. Ship the rotor.

8 Disposal


- i** Electronics are installed in the corresponding rotor body, the remaining accessory parts are free of electrical elements. This rotor body contains only an RFID chip and no batteries or rechargeable batteries.

8.1 Legal regulations

EU countries

In EU member states, electrical and electronic equipment must be disposed of in accordance with Directive 2012/19/EU. This directive has been passed into national law by all EU member states.

Electrical and electronic equipment which has been put on the market after August 13, 2005 must be marked in a special way. According to the European standard the following symbol can be used to mark this equipment:

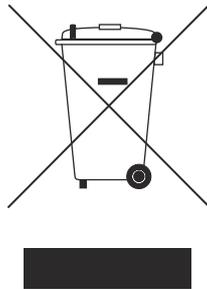

In the EU member states, batteries and rechargeable batteries must be disposed of in accordance with Directive 2006/66/EC. This directive has been passed into national law by all EU member states.

Non-EU countries

Non-EU countries have country-specific standards for the disposal of waste electrical and electronic equipment and the disposal of batteries and rechargeable batteries.

Notes on disposal of electrical and electronic equipment in the United Kingdom

In the United Kingdom, the disposal of electrical and electronic equipment is governed by national regulations which are based on national legislation from 2013, The Waste Electrical and Electronic Equipment Regulations 2013 (as amended), which apply to these devices.

According to these regulations, any electrical and electronic equipment that was put on the market after August 13, 2005 in the business-to-business sector – which applies to this product – must no longer be disposed of with household waste. They are marked with the following symbol to indicate this:



As the disposal regulations may differ from one country to another, please contact your supplier for more information.

8.2 Preparing for disposal

Preparing disposal according to legal regulations

i For information on the legal regulations that apply in your country, contact your local authority and your Eppendorf partner.

1. Check which legal regulations apply to disposal in your country.
2. Choose a certified waste disposal company or contact your Eppendorf partner.
3. Decontaminate and clean the rotor before disposal.
4. If decontamination is not possible, dispose of the rotor as hazardous waste.

Creating a decontamination certificate

Prerequisites:

- The rotor has been decontaminated and cleaned.

1. Download a decontamination certificate from the webpage www.eppendorf.com.
2. Complete the decontamination certificate.

9 Technical data

9.1 Application parameters

Device	Maximum <i>g</i> -force	Maximum speed
SpinPro 6 R	20014 × <i>g</i>	11110 rpm
SpinPro 7 R	25012 × <i>g</i>	12420 rpm

9.2 Service life

i As soon as the maximum number of cycles has been reached, the rotor must be replaced.

Eppendorf SE indicates the maximum service life for rotor and accessories in cycles.

A cycle counts as every centrifugation run in which the rotor is accelerated and then braked again, regardless of speed and duration of the centrifugation run.

The specified service life applies if the following conditions are met:

- Proper use
- Recommended care
- Good condition, free from damage

Components and accessories	Service life from first use
Rotor body	70000 cycles
Aerosol-tight QuickLock Pro rotor lid without replaceable seal	27000 cycles

Components and accessories	Number of autoclaving cycles	Time / temperature
Rotor body	100 autoclaving cycles	20 min / 121 °C
	50 autoclaving cycles	20 min / 134 °C
Aerosol-tight QuickLock Pro rotor lid without replaceable seal	50 autoclaving cycles	20 min / 121 °C

10 Ordering information

Description	Order no.
Rotor FA-10×50 ID Rotor FA-10×50 ID including rotor lid, aerosol-tight 1 pc.	5595 703 002
Rotor lid For rotorFA-10×50 ID, aerosol-tight 1 pc.	5595 704 009

Appendix

Rotor chemical resistance overview

A	Assessment criteria	24
B	Acids and alkalis	25
C	Organic solvents	26
D	Cleaning and decontamination agents	27
E	Other solutions	29

A Assessment criteria

Marking	Resistance	Explanations
■■■	Resistant	Chemical can be used.
■■	Partially resistant	Chemical can cause discoloration of the surface.
■	Not resistant	Chemical can cause severe discoloration or damage to the surface. The chemical must only be used with great caution.
N/A	–	Chemical has not been tested.

B Acids and alkalis

Designation	Concentration	Rotor body and rotor lid	Rotor lid seal	Unlocking mechanism
Ammonia solution	2 %	■ ■	■ ■ ■ ■	■ ■ ■ ■
Ammonia solution	25 %	■	■ ■ ■ ■	■ ■ ■ ■
Acetic acid	12 %	■ ■	■ ■ ■ ■	■ ■ ■ ■
Acetic acid	96 %	■ ■	■ ■ ■ ■	■ ■ ■ ■
Sodium hydroxide	4 %	■	■ ■ ■ ■	■ ■ ■ ■
Sodium hydroxide	20 %	■	■ ■ ■ ■	■ ■ ■ ■
Perchloric acid	10 %	■	■ ■ ■ ■	■
Nitric acid	6.3 %	■	■ ■ ■ ■	■
Nitric acid	65 %	■	■	■
Hydrochloric acid	32 %	■	■ ■ ■ ■	■
Hydrochloric acid	3.6 %	■	■ ■ ■ ■	■ ■ ■ ■
Sulfuric acid	16 %	■	■ ■ ■ ■	■
Sulfuric acid	96 %	■	■	■
Trichloroacetic acid	10 %	■	■ ■ ■ ■	■ ■
Trichloroacetic acid	40 %	■	■ ■ ■ ■	■
Trifluoroacetic acid	10 %	■ ■	■ ■ ■ ■	■ ■

C Organic solvents

Designation	Concentration	Rotor body and rotor lid	Rotor lid seal	Unlocking mechanism
Acetone	≥99.8 %	■■■■	■	■■■■
Acetonitrile	≥99.9 %	■■	■■■■	■■■■
Chloroform	—	■■■■	■■■■	■■■■
Dichloromethane (methylene chloride)	≥99.5 %	■■	■■■■	■■■■
Diethyl ether	—	■■	■■■■	■■
DMSO	10 %	■■	■■■■	■■■■
DMSO	50 %	■■	■■■■	■■■■
DMSO	100 %	■■	■■■■	■■■■
Acetic acid ethyl ester	≥99.5 %	■■■■	■	■■■■
Ethanol	96 %	■■■■	■■■■	■■■■
Formaldehyde	37 %	■■■■	■■■■	■■■■
Isoamyl alcohol	> % 98	■■■■	■■■■	■■■■
Isopropanol	99.8 %	■■■■	■■■■	■■■■
Methanol	99.9 %	■■	■■■■	N/A
Petroleum ether	—	N/A	■■■■	■■■■
Phenol	—	■■	■■■■	■■■■
Toluol	—	■■■■	■■■■	■■■■

D Cleaning and decontamination agents

Designation	Concentration	Rotor body and rotor lid	Rotor lid seal	Unlocking mechanism
Biocidal ZF	—	■■■■	■■■■	■■■■
CaviCide	—	N/A	N/A	N/A
ClearKlens Bi-Spore, VH26	Mix 1 + 48 + 1	N/A	N/A	■■■■
ClearKlens Cleansinald (SC), VH9	0.5 %	N/A	N/A	■■■■
ClearKlens Tego 2000 (SC), VH25	1 %	N/A	N/A	■■■■
COUNT-OFF Liquid Concentrate	2 %	■■■■	■■■■	■■■■
COUNT-OFF Surface Cleaner	—	■■■■	■■■■	■■■■
Descoton Forte	5 %	N/A	N/A	■■■■
Dismozon pur (peroxide-based)	4 %	■■■■	■■■■	■■■■
DNA AWAY	—	■■■■	■■■■	■■■■
DNA-ExitusPlus	—	■■■■	■■■■	■■■■
Ethanol	70 %	N/A	N/A	N/A
Formaldehyde	6 %	■■■■	■■■■	■■■■
Helipur (phenol-based)	6 %	■■■■	■■■■	■■■■
Hexaquart S (QAV-based ⁽¹⁾)	5 %	■■■■	■■■■	■■■■
Hi-TOR Plus (QAV-based ⁽¹⁾)	—	■■■■	■■■■	■■■■
Incidin Plus	3 %	N/A	N/A	N/A
Isopropanol	70 %	N/A	N/A	N/A
Klercide Low Residue Quat Unit Dose Concentrate	01:50	N/A	N/A	■■■■
Klercide Sporicidal Low Residue Peroxide	—	N/A	N/A	■■■■
Korsolex basic (aldehyde-based)	5 %	■■■■	■■■■	■■■■

Designation	Concentration	Rotor body and rotor lid	Rotor lid seal	Unlocking mechanism
Medichem ID Kombi	3 %	N/A	N/A	■■■
Meliseptol (alcohol-based)	—	■■■	■■■	■■■
Sodium hypochlorite	6 %	■■■	■■■	■■■
Sodium hypochlorite	10 %	■	■■■	■■■
Rheosept	—	N/A	N/A	■■■
RNase AWAY	—	■■■	■■■	■■■
RNase-ExitusPlus	—	■■■	■■■	■■■
Sekusept Plus	1.5 %	N/A	N/A	■■■
Spor-Klenz	—	N/A	N/A	■
Hydrogen peroxide	6 %	■■■	■■■	■■■
Hydrogen peroxide	35 %	■■■	■■■	■■■

E Other solutions

Designation	Concentration	Rotor body and rotor lid	Rotor lid seal	Unlocking mechanism
Cesium chloride	1.86g/mL	■■■■	■■■■	■■■■
EDTA (pH 8)	0.5 mol/L	■■■■	■■■■	■■■■
Ficoll (Polysaccharide)	1.077 g/mL	■■■■	■■■■	■■■■
Formamide	50 %	N/A	N/A	N/A
Glutaraldehyde	25 %	■■	■■■■	■
Glycerol	50 %	■■■■	■■■■	■■■■
Guanidinium thiocyanate	4 mol/L	■■	■■■■	■■■■
Guanidinium chloride	6 mol/L	■■■■	■■■■	■■■■
	10 %	N/A	—	—
Mineral oil	—	■■■■	■■■■	■■■■
Sodium acetate (pH 5.2)	2 mol/L	■■■■	■■■■	■■■■
Paraffin oil	—	N/A	N/A	N/A
SDS	1 %	■■	■■■■	■■■■
TRIS buffer pH 7.5	1 mol/L	■■■■	■■■■	■■■■
Triton X-100	1 %	■■■■	■■■■	■■■■
Tween 20	1 %	■■■■	■■■■	■■■■
Water	—	■■■■	■■■■	■■■■

