

Description of changes:

Software version 3.1.6 offers some minor fixes for the Eppendorf ThermoMixer® C, ThermoMixer® F, and ThermoStat® C such as

- incorrect display of characters in the status bar
- translation inconsistency
- rare display glitch during program loading

Moreover, the software version 3.1.6 allows an improved communication between Eppendorf ThermoMixer® C, ThermoMixer® F, and ThermoStat® C and VisioNize® software when connected via the VisioNize box.

Prerequisites:

- Computer with Windows 7 or Windows 10 (Administrator rights are not necessary)
- 5382 000.201 VisioNize Box Adapter cable (or an USB 2.0 high speed A to B device cable shielded)


Instructions:

Determine which software version is currently installed on your device:

- Turn on your Eppendorf ThermoMixer C, ThermoMixer F, or ThermoStat C
- View the software version while the instrument boots

If the installed software is lower than the latest version listed above, perform the upgrade as listed below. You may want to print this page for easy access to the instructions while you work.


1. Select and open the .zip file that was downloaded from eppendorf.com
2. The firmware updates are saved as self-extracting executable files on the desktop
3. Switch on and connect the ThermoMixer C, ThermoMixer F, or ThermoStat C to the USB cable
4. Connect USB cable with the USB port of the computer
5. Double-click on 3.1.6-Eppendorf_ThermoMixer_Autoupdate.exe
6. Allow execution of software update
7. The update starts immediately

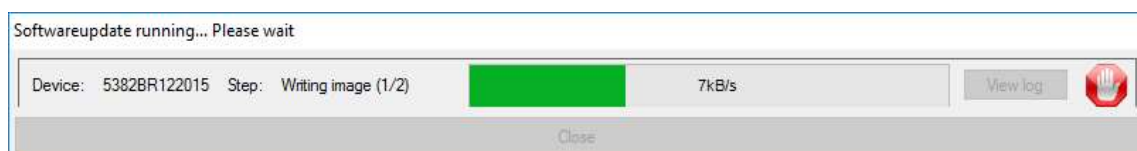

8. Do not interrupt the update process, wait until you are prompted to close the window



9. After the update is completed – you will be asked whether you want to update additional devices

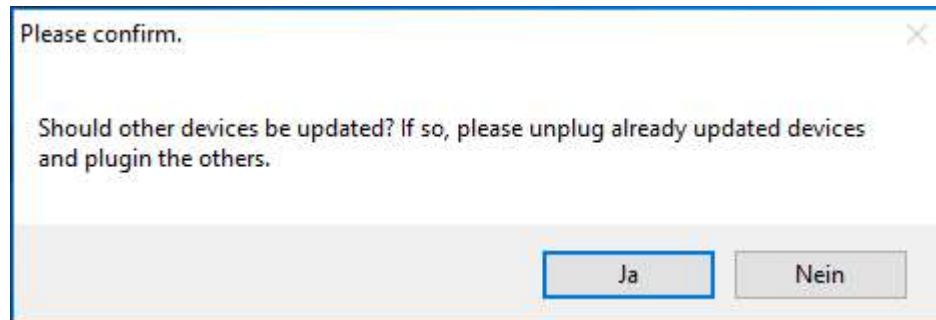

10. Continue with further updates if required

Congratulations! You have installed the latest performance enhancements. Check the Eppendorf website regularly for new revisions.